



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

BS4
9.28.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/22/11 Cal 1 0.103 AC Cal 2 0.101 AC

Between the dates of N/A and 9/22/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on 9/22/11.

Sarah Cione

Sarah L. Cione

9/22/11


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 22, day of September, 20 11.

A.M. Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

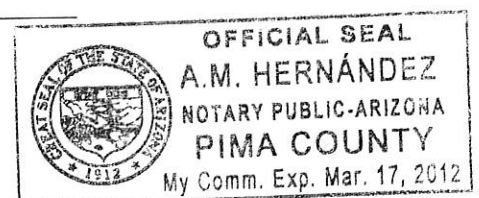

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

B4
9.28.2011

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 9/22/11 TIME 1345

INTOXILYZER SERIAL # 80-000323 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0.100 AC
Result: 0.101 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE

Sarah Cione

2 of 12_{sc}

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

Core Version: 8105.48

09/22/2011

13:47:17

B4

Standard Lot#: 20910100A1

Last Changed By: S. CIONE

#50616

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1980

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:48:21
Diagnostic Test	Pass	13:48:55
Air Blank	0.000	13:49:22
0.100 Cal Check	0.102	13:49:43
Air Blank	0.000	13:50:12
Subject Test	0.000	13:50:37
Air Blank	0.000	13:51:04
Wait	ABT*	
Air Blank	0.000	13:51:41

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 12_{sc}

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

Core Version: 8105.48

09/22/2011

13:52:25

Standard Lot#: 20910100A1

Last Changed By: S. CIONE

#50616

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1980

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:53:12
Diagnostic Test	Pass	13:53:46
Air Blank	RFI*	13:53:56
Air Blank	IPS**	13:54:06

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 12_{sc}

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
09/22/2011

13:54:15

Bh

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:55:30
Diagnostic Test	Pass	13:56:05
Air Blank	0.000	13:56:32
0.100 Cal Check	0.102	13:56:53
Air Blank	0.000	13:57:22
Subject Test	0.047*	14:00:30
Air Blank	0.000	14:01:02
Air Blank	0.000	14:01:29
Subject Test	XXX**	14:01:54
Air Blank	0.000	14:02:27

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5 of 12 sc

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

Core Version: 8105.48

09/22/2011

14:02:51

Standard Lot#: 20910100A1

Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616

TUCSON PD

Operator: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:03:21
0.100 Cal Check	0.101	14:03:42
Air Blank	0.000	14:04:11

Preventative Maintenance
Performed -- Timer Reset

134

EXHIBIT G-2

6 of 12_{sc}

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

154
9.28.2011

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 9/22/11 TIME 1345
Intoxilyzer Serial # 80-000323 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 103 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
09/22/2011

13:45:28

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:45:57
0.100 Cal Check	0.103	13:46:18
Air Blank	0.000	13:46:48

Preventative Maintenance
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

8 of 12 sc

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

BY
9-28-2011

Date: 9/19/11

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-000323

New Lot Number: 20910100A1

New Expiration Date: 8/1/2012

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START".
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change." Select "D" for dry gas and enter the following information:
 - Enter Target Value: 0.100
 - Enter Control Lot #: From the tank label
 - Enter Exp Date: From the tank label (mm/dd/yyyy)
 - Edit info Y/N? Y or N
 - Update Std Info? Y ENTER
 - "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "Control Testing." Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 102 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

9 of 12 sc

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
09/19/2011

14:23:45

dh

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:24:15
0.100 Cal Check	0.102	14:24:36
Air Blank	0.000	14:25:06

Preventative Maintenance
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD

10 of 12 sc

134
9.28.2011

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-000323

Date Returned to Service: 9/22/11



316 E. 9th St., Owensboro KY 42303
Phone 866-835-0690
Fax: 270-685-6268

327183

DATE 8-8-11
RCVD VIA FEDEX

BILL TO:

ATTN:

CITY STATE ZIP

CUSTOMER NO.

PHONE

SHIP TO:

ATTN:

CITY

STATE

ZIP

PHONE

FAX

BILL CODE

MODEL # 2000

SERIAL #

EXTRA PARTS RCVD

DESCRIPTION OF PROBLEM

Infrared / Fuel Cell

BATTERY- EXT N/A V; INT 29 V

VERIFIED PROBLEM?

WORK PERFORMED:

Miscellaneous
PRELIMINARY TEST INFORMATION

WORK PERFORMED:

PARTS USED

[illegible]

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

REPAIR	CAL. ADJ.	CALIBRATION	FINAL
WIV	DR	DR	DZ
1/2 HOURS	HOURS	HOURS	HOURS
3	1	.5	.5
DATE	DATE	DATE	DATE
7/17/11	Aug 30, 2011	Aug 30, 2011	Aug 30, 2011

12 of 12 sc

B4

Certificate of Calibration

Intoxilyzer® 8000

This is to certify that the calibration of **Intoxilyzer®** serial number 80-000323, manufactured by CMI Inc., a subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to be in compliance with the National Highway Traffic Safety Administration Standard for Devices to Measure Breath Alcohol (F.R., No. 179 48705-48710 Sept. 17, 1993) and with International Recommendation OIML R126 Edition 1998 (E). Calibration solutions are traceable to NIST standard material 1828.

Date Aug. 30, 2011

Signed

Dale Rehner

Technician

CMI INC.

316 East 9th Street
Owensboro, KY 42303
USA

Part No. 650519 CMI 7/9/09

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.