



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

154

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/19/07 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 12/22/06 and 1/19/07 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Jennifer L. Coombs
Jennifer L. Coombs

2/9/07

ADPS QAS Permit # 33120

SUBSCRIBED AND SWORN to before me this 9, day of February, 2007.

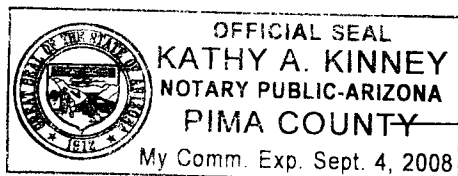

Kathy A. Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Jennifer Coombs 46160 AGENCY Tucson Police DepartmentDATE January 19, 2007 TIME 1516INTOXILYZER SERIAL # 80-001108 LOCATION ODM(☒) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (☒) 1. Clock time check.
(☒) 2. Date check.

OPERATIONAL TESTS

- (☒) 1. Alcohol-free subject test result 0. 000 AC.
(☒) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(☒) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(☒) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

Changed to Mountain Standard Time from
Mountain Daylight
Gas tank change prior to G-2

SIGNATURE

Jennifer L. Coombs

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001108

01/19/2007

15:17:50

Standard Lot#: 629202I

Last Changed By: J. COOMBS

#46160

QAS: J. COOMBS

#46160

TUCSON PD CRIME LAB

Operator: J. COOMBS

#46160

TUCSON PD CRIME LAB

Subject: SQAP,TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:18:59
Diagnostic Test	Pass	15:19:27
Air Blank	0.000	15:19:54
0.100 Cal Check	0.099	15:20:13
Air Blank	0.000	15:20:43
Subject Test	0.000	15:20:58
Air Blank	0.000	15:21:25
Five Minute Wait	ABT*	
Air Blank	0.000	15:22:02

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001108
01/19/2007

15:22:16

Standard Lot#: 629202I
Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160
TUCSON PD CRIME LAB

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:23:44
Diagnostic Test	Pass	15:24:11
Air Blank	0.000	15:24:39
0.100 Cal Check	0.093	15:24:58
Air Blank	0.000	15:25:28
Subject Test	0.000	15:25:51
Air Blank	IPS*	15:25:59
Air Blank	0.000	15:26:27

*Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001108

01/19/2007

15:28:32

Standard Lot#: 629202I

Last Changed By: J. COOMBS

#46160

QAS: J. COOMBS

#46160

TUCSON PD CRIME LAB

Operator: J. COOMBS

#46160

TUCSON PD CRIME LAB

Subject: SQAP,TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:29:29
Diagnostic Test	Pass	15:29:56
Air Blank	0.000	15:30:24
0.100 Cal Check	RFI*	15:30:35
Air Blank	0.000	15:31:06

*RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001108

01/19/2007

15:31:36

Standard Lot#: 629202I

Last Changed By: J. COOMBS

#46160

QAS: J. COOMBS

#46160

TUCSON PD CRIME LAB

Operator: J. COOMBS

#46160

TUCSON PD CRIME LAB

Subject: SQAP,TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:32:29
Diagnostic Test	Pass	15:32:56
Air Blank	0.000	15:33:23
0.100 Cal Check	0.100	15:33:42
Air Blank	0.000	15:34:12
Subject Test	0.270*	15:37:21
Air Blank	0.000	15:37:54
Air Blank	0.000	15:38:22
Subject Test	XXX**	15:38:57
Air Blank	0.000	15:39:31

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001108

01/19/2007

15:41:12


Standard Lot#: 629202I

Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	15:41:40
Cal Check	0.099	15:41:59
Air Blank	0.000	15:42:29

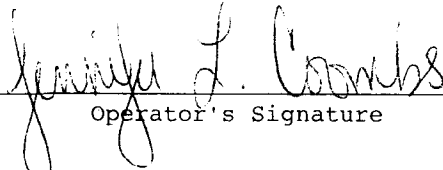

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Jennifer Coombs 46160 AGENCY Tucson Police Department
DATE January 19, 2007 TIME 1516
Intoxilyzer Serial # 80-001108 Location ODM

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Jennifer L. Coombs

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001108

01/19/2007

15:16:25


Standard Lot#: 629202I

Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	15:16:53
Cal Check	0.100	15:17:11
Air Blank	0.000	15:17:42

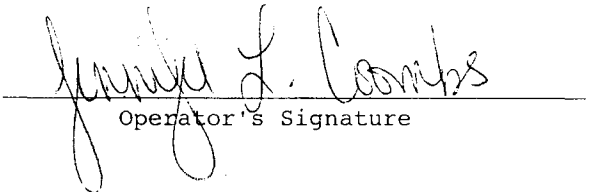
Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

1 of 1
file
B4

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

Date: 1/19/07
QAS: Jennifer Coombs 46160
Intoxilyzer Serial Number: 80-001108
New Lot Number: 629202 I
New Expiration Date: 10/23/2008

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
Enter Target Value: 0.100
Enter Control Lot #: From the tank label
Enter Exp Date: From the tank label (mm/dd/yyyy)
Edit info Y/N? Y or N
Update Std Info? Y ENTER
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE

Jennifer L. Coombs