



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

77


I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3.26.2007 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 2.26.2007 and 3.26.2007 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on                     .  
The instrument was returned to service on                     .





Terry S. Gallegos

April 6, 2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 6<sup>th</sup> day of April, 2007.

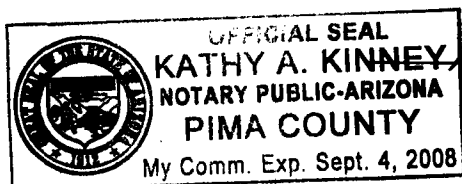

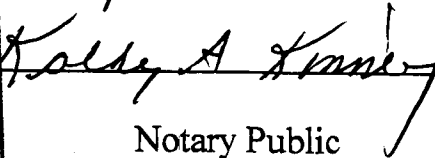  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 8  
TSB  
1/17

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28947 AGENCY Tucson Police Department

DATE 3.26.2007 TIME 1744

INTOXILYZER SERIAL # 80-001163 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.  
( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.  
( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed ✓✓✓  
( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ✓ NO    

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Terry Gallegos

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001163  
03/26/2007

17:46:00

2 of 8  
Blynn

---

Standard Lot#: 635301I  
Last Changed By: G. EPPLEY #26227

---

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	17:46:47
Diagnostic Test	Pass	17:47:14
Air Blank	0.000	17:47:42
0.100 Cal Check	0.099	17:48:00
Air Blank	0.000	17:48:30
Subject Test	0.000	17:48:49
Air Blank	0.000	17:49:16
Five Minute Wait	ABT*	
Air Blank	0.000	17:49:52

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

3 of 8  
TS6mm

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001163  
03/26/2007

17:49:56

---

Standard Lot#: 635301I  
Last Changed By: G. EPPLEY #26227

---

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	17:50:44
Diagnostic Test	Pass	17:51:12
Air Blank	RFI*	17:51:30
Air Blank	0.000	17:51:58

\*RFI Detect

### **Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001163  
03/26/2007

17:52:18

4 of 8  
By [signature]

---

Standard Lot#: 635301I  
Last Changed By: G. EPPLEY #26227

---

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	IPS*	17:53:10
Air Blank	0.000	17:53:39

\*Improper Sample

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001163  
03/26/2007

17:53:58

50f8  
Dy 11/11

Standard Lot#: 635301I  
Last Changed By: G. EPPLEY #26227

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:54:49
Diagnostic Test	Pass	17:55:16
Air Blank	0.000	17:55:43
0.100 Cal Check	0.099	17:56:02
Air Blank	0.000	17:56:32
Subject Test	0.097*	17:59:40
Air Blank	0.000	18:00:12
Air Blank	0.000	18:00:39
Subject Test	XXX**	18:00:54
Air Blank	0.000	18:01:27

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001163

03/26/2007

18:01:54

60f 8  
TS4 mm

-----  
Standard Lot#: 635301I

Last Changed By: G. EPPLEY #26227  
-----

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	18:02:21
Cal Check	0.099	18:02:40
Air Blank	0.000	18:03:10

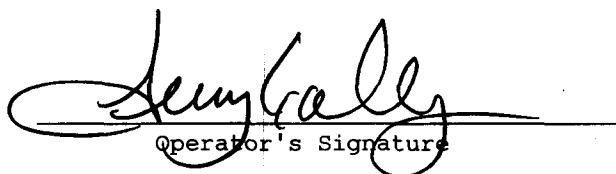
  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos AGENCY Tucson Police Department  
DATE 3.26.2007 TIME 1744  
Intoxilyzer Serial # 80-001163 Location DAI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)



INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001163

03/26/2007

17:44:35

-----  
Standard Lot#: 635301I

Last Changed By: G. EPPLEY #26227  
-----

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	17:45:02
Cal Check	0.099	17:45:21
Air Blank	0.000	17:45:50

  
Operator's Signature

80f8  
By [signature]



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

10f2  
TS4  
mm

Date: 3/22/07  
QAS: OFC. EAPLEY 26227  
Intoxilyzer Serial Number: 80-001163  
New Lot Number: 635301 I  
New Expiration Date: 122808

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
Enter Target Value: 0.100  
Enter Control Lot #: From the tank label  
Enter Exp Date: From the tank label (mm/dd/yyyy)  
Edit info Y/N? Y or N  
Update Std Info? Y ENTER  
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE

[Signature] # 26227


2 of 2  
T84  
NW

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-081163  
03/22/2007 03:23:35

Standard Lot#: 6353011  
Last Changed By: G. EPPLEY #26227

DAS: G. EPPLEY #26227  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	03:24:03
Cal Check	0.100	03:24:21
Air Blank	0.000	03:24:51

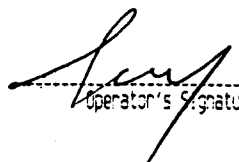  
Operator's Signature