



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001199 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
11-7-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/14/11 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 9/15/11 and 10/14/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


10/27/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 27, day of October, 20 11.

A.M. Hernández

Notary Public

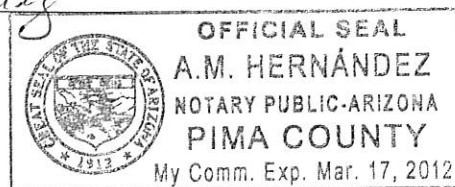

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

1 of 7_{sc}

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

By
11.7.2011

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 10/14/11 TIME 1007

INTOXILYZER SERIAL # 80-001199 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001199

10/14/2011

10:10:12

B4

Standard Lot#: 20910100A1

Last Changed By: S. CIONE

#50616

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:10:58
Diagnostic Test	Pass	10:11:26
Air Blank	0.000	10:11:53
0.100 Cal Check	0.099	10:12:12
Air Blank	0.000	10:12:41
Subject Test	0.000	10:13:09
Air Blank	0.000	10:13:36
Five Minute Wait	ABT*	
Air Blank	0.000	10:14:24

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001199

10/14/2011

10:14:56

TS4

Standard Lot#: 20910100A1

Last Changed By: S. CIONE

#50616

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:15:27
Air Blank	IPS**	10:15:35

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7_{sc}

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001199

10/14/2011

10:17:28

T84

Standard Lot#: 20910100A1

Last Changed By: S. CIONE

#50616

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:18:31
Diagnostic Test	Pass	10:18:58
Air Blank	0.000	10:19:26
0.100 Cal Check	0.100	10:19:44
Air Blank	0.000	10:20:14
Subject Test	0.126*	10:23:22
Air Blank	0.000	10:23:54
Air Blank	0.000	10:24:21
Subject Test	XXX**	10:25:07
Air Blank	0.000	10:25:40

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001199

10/14/2011

10:26:28


Standard Lot#: 20910100A1

Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616TUCSON PD

TS4

Test	g/210L	Time
Air Blank	0.000	10:26:56
Cal Check	0.098	10:27:14
Air Blank	0.000	10:27:44

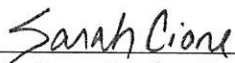

Operator's Signature

6 of 7sc

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

T34
11.7.2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 10/14/11 TIME 1007
Intoxilyzer Serial # 80-001199 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7_{sc}

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001199
10/14/2011 10:07:40

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

TS4

Test	g/210L	Time
Air Blank	0.000	10:08:07
Cal Check	0.099	10:08:26
Air Blank	0.000	10:08:55

Sarah Cione
Operator's Signature