



# CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

B4  
6-7-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002262 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/26/10 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 4/28/10 and 5/26/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione

6/3/10


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 3rd, day of June, 2010.

A.M. Hernández

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

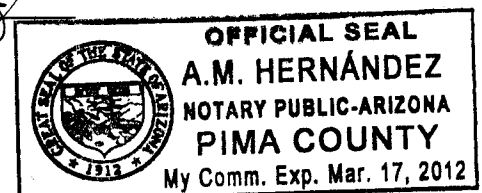

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4  
6-7-2010

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 5/26/10 TIME 857

INTOXILYZER SERIAL # 80-002262 LOCATION Team 2 - DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0.000 AC.

(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0.100 AC  
Result: 0.099 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Sarah Cione

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
05/26/2010

08:59:43

2 of 7<sub>sc</sub>

dy

---

Standard Lot#: 642845  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:00:29
Diagnostic Test	Pass	09:00:56
Air Blank	0.000	09:01:24
0.100 Cal Check	0.100	09:01:42
Air Blank	0.000	09:02:12
Subject Test	0.000	09:02:53
Air Blank	0.000	09:03:21
Five Minute Wait	ABT*	
Air Blank	0.000	09:04:02

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7<sub>sc</sub>

B4

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
05/26/2010

09:04:44

---

Standard Lot#: 642845  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	09:05:08
Air Blank	IPS**	09:05:16

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002262  
05/26/2010

09:05:22

4 of 7<sub>sc</sub>

154

---

Standard Lot#: 642845  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:06:23
Diagnostic Test	Pass	09:06:51
Air Blank	0.000	09:07:18
0.100 Cal Check	0.100	09:07:37
Air Blank	0.000	09:08:06
Subject Test	0.036*	09:11:14
Air Blank	0.000	09:11:46
Air Blank	0.000	09:12:13
Subject Test	XXX**	09:12:38
Air Blank	0.000	09:13:11

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002262

05/26/2010

09:14:07

-----  
Standard Lot#: 642845

Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:14:35
Cal Check	0.099	09:14:53
Air Blank	0.000	09:15:23

*Sarah Cione*

Operator's Signature

5 of 7<sub>sc</sub>

TS4

## EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 75c

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000TS4  
6-7-2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 5/26/10 TIME 857  
Intoxilyzer Serial # 80-002262 Location Team 2 - DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002262  
05/26/2010

08:57:02

-----  
Standard Lot#: 642845  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:57:29
Cal Check	0.100	08:57:48
Air Blank	0.000	08:58:17

Sarah Cione  
Operator's Signature

7 of 7<sub>sc</sub>

154