



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002803 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSY
5.2.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4/28/11 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 3/31/11 and 4/28/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -

The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


4/28/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 28th day of April, 20 11.

A.M. Hernández

Notary Public

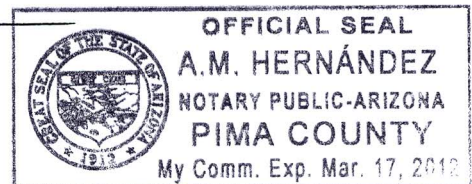

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7_{sc}

154
5.2.2011

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 4/28/11 TIME 1057

INTOXILYZER SERIAL # 80-002803 LOCATION Team 2

(☒) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (☒) 1. Clock time check.
(☒) 2. Date check.

OPERATIONAL TESTS

- (☒) 1. Alcohol-free subject test result 0. 000 AC.
(☒) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(☒) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(☒) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002803
04/28/2011

10:59:09

2 of 7sc
Bh

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:59:50
Diagnostic Test	Pass	11:00:18
Air Blank	0.000	11:00:45
0.100 Cal Check	0.097	11:01:04
Air Blank	0.000	11:01:33
Subject Test	0.000	11:01:52
Air Blank	0.000	11:02:19
Five Minute Wait	ABT*	
Air Blank	0.000	11:02:59

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002803
04/28/2011

11:03:06

3 of 7_{sc}

154

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	11:03:28
Air Blank	IPS**	11:03:37

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002803
04/28/2011

11:03:42

4 of 7_{sc}
BH

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:04:23
Diagnostic Test	Pass	11:04:51
Air Blank	0.000	11:05:18
0.100 Cal Check	0.097	11:05:37
Air Blank	0.000	11:06:07
Subject Test	0.050*	11:09:15
Air Blank	0.000	11:09:48
Air Blank	0.000	11:10:15
Subject Test	XXX**	11:10:32
Air Blank	0.000	11:11:05

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002803
04/28/2011 11:11:20

5 of 7sc

154

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:11:47
Cal Check	0.097	11:12:06
Air Blank	0.000	11:12:36



Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

by
5.2.2011

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 4/28/11 TIME 1057
Intoxilyzer Serial # 80-002803 Location Team 2

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002803

04/28/2011

10:57:49

Standard Lot#: 719919

Last Changed By: S. CIONE #50616

OAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:58:16
Cal Check	0.097	10:58:35
Air Blank	0.000	10:59:05

Sarah C. Cione
Operator's Signature

7 of 7
sc

154