



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 5000 Serial Number 64-002516 (#12) maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R9-14-404.A.3, for a 0.080AC solution are 0.070AC to 0.090AC.

Date: 2.15.06 Cal 1 0.078 AC Cal 2 0.080 AC

Date: 3.16.06 Cal 1 0.075 AC Cal 2 0.074 AC

Date: _____ Cal 1 _____ AC Cal 2 _____ AC

Between the dates of 2.15.06 and 3.16.06, the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

Terry S. Gallegos


Terry S. Gallegos

3.23.2006

ADPS QAS Permit # QAS 976

SUBSCRIBED AND SWORN to before me this 23, day of

March, 2006.

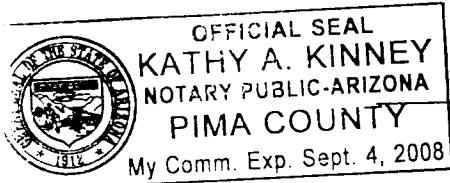

Kathy A. Kinney

Notary Public

*Includes documents specified in A.A.C. R9-14-404.A.6

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT Q

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404(A)

10f6
TS4

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

AGENCY Tucson Police QA SPECIALIST Terry Gallegos 28847
 (Printed Name)
 INTOXILYZER SERIAL # 64-002516 (#12) LOCATION Traffic
 DATE March 16 TIME 1904

✓ Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

DIAGNOSTIC TESTS

- DVM Test check. Setting should be between .200 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is .334.
- Display Test check. Mode selection switch S1 on and S2 and S3 off or keyboard menu selection "V".
- Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "P".
- Clock time check. Mode selection switch S10 on or keyboard menu selection "E".
- * Date check. Mode selection switch S11 on or keyboard menu selection "E".

OPERATIONAL TESTS

- Alcohol free subject Test result 0.000 AC
- Error recognition logic system functioning
 - Invalid test printed
- Proper sample recognition system
 - Invalid sample printed
 - Deficient sample printed
- N/A Completeness of sample purge with collection tube**
 - Calibration standard 0.080 AC
 - Results: 0.075AC and 0.074 AC

Instrument operating properly and accurately. YES NO

COMMENTS: Changed date from 3/15/06 to 3/16/06

SIGNATURE

0. N/A alcohol concentration standard collected for subsequent analysis.**

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

**THIS STEP IS ONLY REQUIRED IF THIS DEVICE IS BEING USED FOR
SAMPLE CAPTURE.

#12		
	SUBJECT'S NAME Traffic	INSTRUMENT LOCATION
		OPERATOR
		ADDITIONAL INFORMATION AND / OR REMARKS
		Gallegos 28847 SERVED

SUBJECT'S NAME	Traffic	
INSTRUMENT LOCATION		
TIME FIRST OBSERVED	Gallegos 28847	
OPERATOR		
ADDITIONAL INFORMATION AND/OR REMARKS		

SUBJECT'S NAME	Big Traffic	
TIME FIRST OBSERVED		
INSTRUMENT LOCATION		
ON TRAILOR	Galleys 28847	
ADDITIONAL INFORMATION AND/or REMARKS		

TEST	BAC	VALUE	TIME
AIR BLANK	.000		19:04
SUBJECT TEST	.000		19:05
AIR BLANK	.000		19:05

GMI INC

EMI INC.

COMMERCIAL INSTRUMENT PRINTER CARD

INTOXILYZER® INSTRUMENT PRINTER CARD

INTOXICATOR

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

Tucson Police Department Crime Laboratory

FEDERAL SIGNAL CORP. / CMI INC
 INTOXILYZER - ALCOHOL ANALYZER
 AZ MODEL 5000 SN 64-002516
 03/15/86

TEST	BAC VALUE	TIME	TEST	BAC VALUE	TIME
AIR BLANK	.000	19:11	AIR BLANK	.000	19:13
CAL. CHECK	.075	19:11	CAL. CHECK	.074	19:13
AIR BLANK	.000	19:12	AIR BLANK	.000	19:13

#12

SUBJECT'S NAME	<u>Gallenos</u>
INSTRUMENT LOCATION	<u>Traffic</u>
TIME FIRST OBSERVED	<u>28847</u>
OPERATOR	
ADDITIONAL INFORMATION AND/OR REMARKS	

30feb86
JL

SUBJECT'S NAME	<u>Gallenos</u>
INSTRUMENT LOCATION	<u>Traffic</u>
TIME FIRST OBSERVED	<u>28847</u>
OPERATOR	
ADDITIONAL INFORMATION AND/OR REMARKS	

CMI INC

INTOXILYZER INSTRUMENT PRINTER CARD
 © 1986 by CMI INC

CMI INC

INTOXILYZER INSTRUMENT PRINTER CARD
 © 1986 by CMI INC

404
TSL

SN 64-002516 03/15/86
E47.02 19:10
INVALID TEST

#12

SUBJECT'S NAME		
Traffic		
TIME FIRST OBSERVED	INSTRUMENT LOCATION	
Gallegos 28847		
OPERATOR		
ADDITIONAL INFORMATION AND/OR REMARKS		
Pushed Button wrong time		
INTOXILYZER® INSTRUMENT PRINTER CARD		
© 1986 by CMI INC.		

CMI

EXHIBIT P

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404(A)

*50f6
B4*

ARIZONA DEPARTMENT OF HEALTH SERVICES


STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

STANDARD CALIBRATION CHECK PROCEDURE

Agency Tucson Police Department Date 3.16.2006 Time 1904
 Intoxilyzer Serial # 64-002516 (#12) Location Traffic
 QA Specialist Terry Gallegos 28847
 (Print Name)

- (✓) 1. Pour a standard alcohol solution of known value into a clean dry simulator jar and assemble the simulator. Insure that a tight seal has been made.
Standard value: 0.080 AC
- (✓) 2. Turn on the simulator and allow the temperature to reach $34^{\circ}\text{C} \pm .2^{\circ}\text{C}$.
- (✓) 3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- (✓) 4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- (✓✓) 5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- (✓✓) 6. Push Start Test button or press enter on keyboard.
- (✓✓) 7. Insert card in response to display.
- (✓✓) 8. Air blank completed.
- (✓✓) 9. Calibration check completed.
Test results 0.075AC and 0.074 AC
- (✓✓) 10. Air blank completed.
- (✓✓) 11. When display reads Test Complete remove evidence card. Attach the card to the completed checklist.
- (✓) 12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE

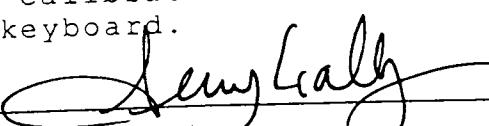

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DHS/DSLS/Form C116 (Rev. 7-93)

FEDERAL SIGNAL CORP. / CMI INC
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SH 64-002516
03/16/86

TEST	BAC	VALUE	TIME
AIR BLANK	.000	19:16	
INVALID SAMPLE	.XXX	19:16	
AIR BLANK	.000	19:16	

FEDERAL SIGNAL CORP. / CMI INC
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SH 64-002516
03/16/86

TEST	BAC	VALUE	TIME
AIR BLANK	.000	19:18	
*SUBJECT TEST	.056	19:20	
AIR BLANK	.000	19:21	

* DEFICIENT SAMPLE - VALUE
PRINTED WAS HIGHEST OBTAINED.

SUBJECT'S NAME		INSTRUMENT LOCATION	
#12		Traffic	
TIME FIRST OBSERVED		OPERATOR	
		ADDITIONAL INFORMATION AND/OR REMARKS	
Gallagos 20847		cmi	
Shut breath		INTOXILYZER™ INSTRUMENT PRINTER CARD	
		© 1986 by CMI INC.	
		INTOXILYZER™ INSTRUMENT PRINTER CARD	
		© 1986 by CMI INC.	

60th
Tdy