



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 5000 Serial Number 64-003336 (#15) maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program. SM

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.080AC solution are 0.070AC to 0.090AC.

Date: 5/17/06 Cal 1 0.079 AC Cal 2 0.081 AC

Date: 6/15/06 Cal 1 0.080 AC Cal 2 0.080 AC

Date: N/A Cal 1 AC Cal 2 AC

Between the dates of 5/17/06 and 6/15/06 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Jennifer L. Coombs
Jennifer L. Coombs

6/20/06

ADPS QAS Permit # QAS 29886

SUBSCRIBED AND SWORN to before me this 20, day of June, 2006.

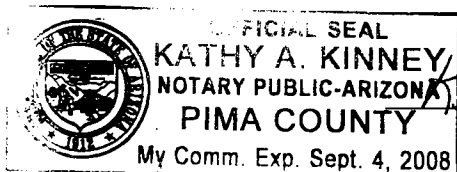

Kathy A. Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Jennifer Coombs 4/6/00 AGENCY Tucson Police Department
DATE June 15, 2006 TIME 1010
INTOXILYZER SERIAL # 64-003336 (#15) LOCATION DUI HQ

- (✓) 1. Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

DIAGNOSTIC TESTS

- (✓) 1. DVM Test check. Setting should be between .010 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is
- .363
- .
-
- (✓) 2. Display Test check. Mode selection switch S1 on and S2 and S3 off or keyboard menu selection "V".
-
- (✓) 3. Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "P".
-
- (✓) 4. Clock time check. Mode selection switch S10 on or keyboard menu selection "E".
-
- (✓) 5. Date check. Mode selection switch S11 on or keyboard menu selection "E".

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject Test result 0.
- 000
- AC
-
- (✓) 2. Error recognition logic system functioning
-
- Invalid test printed
-
- (✓) 3. Proper sample recognition system
-
- Invalid sample printed ✓
-
- Deficient sample printed ✓
-
- (✓) 4. Standard Calibration Check standard 0.
- 080
- AC
-
- Results: 0.
- 080 AC
- & 0.
- 080
- AC

Instrument operating accurately and properly. YES ☒ NO ☐

COMMENTS:

SIGNATURE

Jennifer L Coombs

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DPS Form Exh E-4 (Rev. 05-1)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER™ INSTRUMENT PRINTER CARD
© 1986 by CMI INC.

NTQXIIYZER[®] INSTRUMENT PRINTER CARD

Tucson Police Department Crime Laboratory

3 of 6
file

#15

CMI INC
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 64-003336
06/15/06

TEST	BAC VALUE	TIME
AIR BLANK	.000	10:43
CAL. CHECK	.080	10:43
AIR BLANK	.000	10:43

SUBJECT'S NAME

DUI #10

INSTRUMENT LOCATION

TIME FIRST OBSERVED

46160

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

#15

CMI INC
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 64-003336
06/15/06

TEST	BAC VALUE	TIME
AIR BLANK	.000	10:41
CAL. CHECK	.080	10:42
AIR BLANK	.000	10:42

SUBJECT'S NAME

DUI #10

INSTRUMENT LOCATION

TIME FIRST OBSERVED

46160

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

4 of 6
TS4

#15

SN 64-003336
E47.04
INVALID TEST
SAMPLE INTRODUCED
AT IMPROPER TIME.

06/15/06
10:21

SUBJECT'S NAME

TIME FIRST OBSERVED

INSTRUMENT LOCATION

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

INTOXILYZER[®] INSTRUMENT PRINTER CARD

© 1986 by CMI INC

CMI INC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*STANDARD CALIBRATION CHECK PROCEDURE *bm**5d6*
je

QA SPECIALIST Jennifer Combs 46160 AGENCY Tucson Police Department
DATE June 15, 2006 TIME 1010
INTOXILYZER SERIAL # 64-003336 (#15) LOCATION DUI HQ

- (✓) 1. Pour a standard alcohol concentration solution into a clean dry simulator jar and assemble the simulator. Ensure that a tight seal has been made. Standard value: 0. 080 AC.
- (✓) 2. Turn on the simulator and allow the temperature to reach 34°C ± 0.2°C.
- (✓) 3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- (✓) 4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- ✓ (✓) 5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- ✓ (✓) 6. Push Start Test button or press enter on keyboard.
- ✓ (✓) 7. Insert card in response to display.
- ✓ (✓) 8. Air blank completed.
- ✓ (✓) 9. Standard Calibration Check completed. Test results 0. 080 AC & 0. 080 AC.
- ✓ (✓) 10. Air blank completed.
- ✓ (✓) 11. When test complete, remove printed record. Attach the record to the completed checklist.
- (✓) 12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE

Jennifer L. Combs

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DPS Form Exh E-2 (Rev. 05-1)

#15

CMI INC
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 64-003336
06/15/06

TEST	BAC VALUE	TIME
AIR BLANK	.000	10:22
INVALID SAMPLE	.XXX	10:23
AIR BLANK	.000	10:23

SUBJECT'S NAME

DUL H0

INSTRUMENT LOCATION

TIME FIRST OBSERVED

461600

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

Sweet Breath

INTOXILYZER* INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

CMI INC

#15

CMI INC
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 64-003336
06/15/06

TEST	BAC VALUE	TIME
AIR BLANK	.000	10:24
*SUBJECT TEST	.136	10:26
AIR BLANK	.000	10:26

* DEFICIENT SAMPLE - VALUE
PRINTED WAS HIGHEST OBTAINED.

SUBJECT'S NAME

DUL H0

INSTRUMENT LOCATION

TIME FIRST OBSERVED

461600

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

INTOXILYZER* INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

CMI INC