



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 5000 Serial Number 66-004798 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R9-14-404.A.3, for a 0.080AC solution are 0.070AC to 0.090AC.

Date: 3.16.06 Cal 1 0.078 AC Cal 2 0.079 AC

Date: 4.13.06 Cal 1 0.078 AC Cal 2 0.080 AC

Date: \_\_\_\_\_ Cal 1 \_\_\_\_\_ AC Cal 2 \_\_\_\_\_ AC

Between the dates of 3.16.06 and 4.13.06 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on \_\_\_\_\_.

The instrument was returned to service on \_\_\_\_\_.

Terry S. Gallegos


Terry S. Gallegos

4.19.2006

ADPS QAS Permit # QAS 976

SUBSCRIBED AND SWORN to before me this 19, day of

April, 2006.

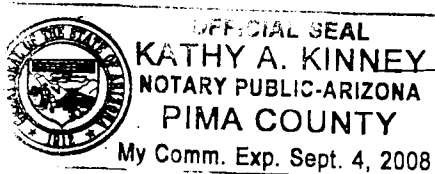

Kathy A. Kinney

Notary Public

\*Includes documents specified in A.A.C. R9-14-404.A.6

## EXHIBIT Q

1076 TS4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 5000\*

AGENCY Tucson Police QA SPECIALIST Terry S. Gallegos 28847  
(Print Name)

INTOXILYZER SERIAL # 66-004798 LOCATION OSD  
DATE April 13 2006 TIME 0907

☒ Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

## DIAGNOSTIC TESTS

☒ DVM Test check. Setting should be between .200 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is .330.

☒ Display Test check. Mode selection switch S1 on and S2 and S3 off or keyboard menu selection "V".

☒ Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "P".

☒ Clock time check. Mode selection switch S10 on or keyboard menu selection "E".

☒ Date check. Mode selection switch S11 on or keyboard menu selection "E".

## OPERATIONAL TESTS

☒ Alcohol free subject Test result 0.000 AC

☒ Error recognition logic system functioning  
Invalid test printed

☒ Proper sample recognition system  
Invalid sample printed ☒  
Deficient sample printed ☒

N/A Completeness of sample purge with collection tube\*\*


☒ Calibration standard 0.080 AC

Results: 0.078 AC and 0.080 AC

Instrument operating properly and accurately. YES ☒ NO ☐

COMMENTS:

SIGNATURE

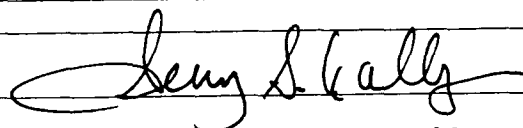

0. N/A alcohol concentration standard collected for subsequent analysis.\*\*

\*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD  
\*\*THIS STEP IS ONLY REQUIRED IF THIS DEVICE IS BEING USED FOR  
SAMPLE CAPTURE.

DHS/DSLS/Form C117 (Rev. 7-93)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

TEST	SRAC VALUE	TIME
AIR BLANK	.000	99:07
SUBJECT TEST	.000	99:30
AIR BLANK	.000	99:00

$$DVM = .330$$

SUBJECT'S NAME \_\_\_\_\_

INSTRUMENT LOCATION \_\_\_\_\_

0DD

Galleos 28847

ADDITIONAL INFORMATION AND/OR REMARKS

ANTOXILYZER<sup>®</sup> INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

[illegible]

SUBJECT'S NAME \_\_\_\_\_

INSTRUMENT LOCATION \_\_\_\_\_

0DD

Gallegos	28847
OPERATOR	
ASTUTICAL INFORMATION AND /OR REMARKS	

ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

20th TS4

# Tucson Police Department Crime Laboratory

3046  
TSH

TUCSON POLICE DEPARTMENT  
INTOXILYZER - ALCOHOL ANALYZER  
AZ MODEL 5000 SN 66-004799  
24/13/96 09:12

TEST  
AIR BLANK 09:11  
CAL. CHECK 09:11  
AIR BLANK 09:11

BRAC VALUE  
0.000  
0.000  
0.000

TUCSON POLICE DEPARTMENT  
INTOXILYZER - ALCOHOL ANALYZER  
AZ MODEL 5000 SN 66-004799  
24/13/96 09:10

TEST  
AIR BLANK 09:09  
CAL. CHECK 09:10  
AIR BLANK 09:10

BRAC VALUE  
0.000  
0.070  
0.000

SUBJECT'S NAME  
0DD

TIME FIRST OBSERVED  
INSTRUMENT LOCATION

Gallegos 28847

OPERATOR  
ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

SUBJECT'S NAME  
0DD

TIME FIRST OBSERVED  
INSTRUMENT LOCATION

Gallegos 28847

OPERATOR  
ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

4056  
B4

je

SN 00-00-180  
EP07.14  
INVALID TEST

04/19/86  
09:23

PRESSED START TEST AT  
IMPROPER TIME

SUBJECT'S NAME	
TIME FIRST OBSERVED	INSTRUMENT LOCATION
Gallagos 20847	
OPERATOR	
ADDITIONAL INFORMATION AND/OR REMARKS	

INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC

**CMI**

## EXHIBIT P

506  
784THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R9-14-404(A)

## ARIZONA DEPARTMENT OF HEALTH SERVICES


STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 5000\*

## STANDARD CALIBRATION CHECK PROCEDURE

Agency Tucson Police Department Date 4.13.2006 Time 0907  
Intoxilyzer Serial # 66-004798 Location ODD  
QA Specialist Terry S. Gallegos 28847  
(Print Name)

- (✓) 1. Pour a standard alcohol solution of known value into a clean dry simulator jar and assemble the simulator. Insure that a tight seal has been made. Standard value: 0.080 AC
- (✓) 2. Turn on the simulator and allow the temperature to reach  $34^{\circ}\text{C} \pm .2^{\circ}\text{C}$ .
- (✓) 3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- (✓) 4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- (✓✓) 5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- (✓✓) 6. Push Start Test button or press enter on keyboard.
- (✓✓) 7. Insert card in response to display.
- (✓✓) 8. Air blank completed.
- (✓✓) 9. Calibration check completed. Test results 0.078AC and 0.080 AC
- (✓✓) 10. Air blank completed.
- (✓✓) 11. When display reads Test Complete remove evidence card. Attach the card to the completed checklist.
- (✓) 12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE

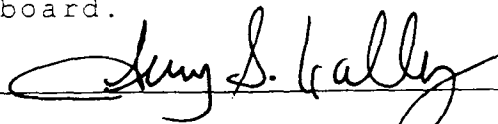

\*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DHS/DSLS/Form C116 (Rev. 7-93)

TUCSON POLICE DEPARTMENT  
INTOXILYZER - ALCOHOL ANALYZER  
AZ MODEL 5000 SH 66-004798  
04/13/86 09:17

TEST	BRAC VALUE	TIME
AIR BLANK	.000	09:17
INVALID SAMPLE	.004	09:17
AIR BLANK	.000	09:17

BEGIN NEW DEPRIVATION PERIOD

SUBJECT'S NAME  
TIME FIRST OBSERVED  
INSTRUMENT LOCATION

Gallegos 29847

OPERATOR  
ADDITIONAL INFORMATION AND/OR REMARKS

Sweet breath



INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

TUCSON POLICE DEPARTMENT  
INTOXILYZER - ALCOHOL ANALYZER  
AZ MODEL 5000 SH 66-004798  
04/13/86 09:21

TEST	BRAC VALUE	TIME
AIR BLANK	.000	09:18
*SUBJECT TEST	.085	09:21
AIR BLANK	.090	09:21

\* DEFICIENT SAMPLE - VALUE  
PRINTED WAS HIGHEST OBTAINED.

SUBJECT'S NAME  
TIME FIRST OBSERVED  
INSTRUMENT LOCATION

Gallegos 29847

OPERATOR  
ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

60f6  
TS4  
fle