



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

154

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 5000 Serial Number 66-005004 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R9-14-404.A.3, for a 0.080AC solution are 0.070AC to 0.090AC.

Date: 11/2/05 Cal 1 0.081 AC Cal 2 0.080 AC

Date: 12/2/05 Cal 1 0.080 AC Cal 2 0.079 AC

Date: N/A Cal 1                      AC Cal 2                      AC

Between the dates of 11/2/05 and 12/2/05 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on 12/2/05.  
The instrument was returned to service on N/A.

Jennifer L. Coombs  
Jennifer L. Coombs

12/2/05

ADPS QAS Permit # QAS 29886

SUBSCRIBED AND SWORN to before me this 2, day of December, 2005.

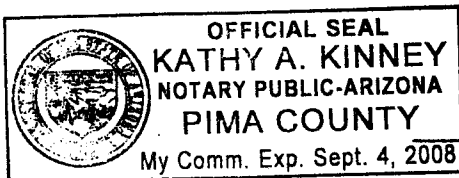

Kathy A Kinney  
Notary Public

\*Includes documents specified in A.A.C. R9-14-404.A.6

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

## EXHIBIT Q

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R9-14-404(A)

## ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES *ph*  
INTOXILYZER MODEL 5000\*AGENCY Tucson Police QA SPECIALIST Jennifer Coombs 46160  
(Print Name)INTOXILYZER SERIAL # 66-005004 LOCATION LAB  
DATE December 2<sup>nd</sup> 2005 TIME 1330

☒ Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

## DIAGNOSTIC TESTS

- ☒ DVM Test check. Setting should be between .200 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is .407.
- ☒ Display Test check. Mode selection switch S1 on and S2 and S3 off or keyboard menu selection "V".
- ☒ Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "P".
- ☒ Clock time check. Mode selection switch S10 on or keyboard menu selection "E".
- ☒ Date check. Mode selection switch S11 on or keyboard menu selection "E".

## OPERATIONAL TESTS

- ☒ Alcohol free subject Test result 0.000 AC
- ☒ Error recognition logic system functioning  
Invalid test printed
- ☒ Proper sample recognition system  
Invalid sample printed ✓  
Deficient sample printed ✓
- N/A Completeness of sample purge with collection tube\*\*
- ☒ Calibration standard 0.080 AC  
Results: 0.080 AC & 0.079 AC
- Instrument operating properly and accurately. YES ☒ NO ☐

COMMENTS: \*moved clock 1336 → 1341

SIGNATURE

Jennifer L. Coombs  
0. N/A alcohol concentration standard collected for subsequent analysis.\*\*

\*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD  
\*\*THIS STEP IS ONLY REQUIRED IF THIS DEVICE IS BEING USED FOR  
SAMPLE CAPTURE.

DHS/DSLS/Form C117 (Rev. 7-93)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

1 of 6  
file

| TEST         | BRAC VALUE | TIME  |
|--------------|------------|-------|
| AIR BLANK    | .000       | 13:47 |
| SUBJECT TEST | .000       | 13:48 |
| AIR BLANK    | .000       | 13:48 |

$$\underline{DVM = .407}$$

**SUBJECT'S NAME**

LAB

| INSTRUMENT LOCATION |
|---------------------|
| LTB                 |

TIME FIRST OBSERVED

46100

OPERATOR

**ADDITIONAL INFORMATION AND/OR REMARKS**



**INTOXILYZER® INSTRUMENT PRINTER CARD**

© 1986 by CMI INC.

[illegible]

SUBJECT'S NAME

AL

INSTRUMENT LOCATION

TIME FIRST OBSERVED

440160

5

ADDITIONAL INFORMATION AND/OR REMARKS:



**INTOXILYZER<sup>®</sup> INSTRUMENT PRINTER CARD**

© 1986 by CMI INC.

3 of 6  
file

# Tucson Police Department Crime Laboratory

|  |  |  |  |                                    |  |                                       |  |
|--|--|--|--|------------------------------------|--|---------------------------------------|--|
| TUCSON POLICE DEPARTMENT<br>INTOXILYZER - ALCOHOL ANALYZER<br>AZ MODEL 5000 SN 66-005004<br>12/02/05 13:47 |  | TEST<br>AIR BLANK<br>CAL. CHECK<br>AIR BLANK |  | BRAC VALUE<br>.000<br>.079<br>.000 |  | TIME<br>13:46<br>13:46<br>13:47       |  |
| TIME FIRST OBSERVED  |  | SUBJECT'S NAME                               |  | LAB                                |  | INSTRUMENT LOCATION                   |  |
| 461160   |  | OPERATOR                                     |  | 461160                             |  | ADDITIONAL INFORMATION AND/OR REMARKS |  |
| INTOXILYZER® INSTRUMENT PRINTER CARD   |  | © 1986 by CMI INC.                           |  | CMI                                |  |                                       |  |

|  |  |  |  |                                    |  |                                       |  |
|--|--|--|--|------------------------------------|--|---------------------------------------|--|
| TUCSON POLICE DEPARTMENT<br>INTOXILYZER - ALCOHOL ANALYZER<br>AZ MODEL 5000 SN 66-005004<br>12/02/05 13:45 |  | TEST<br>AIR BLANK<br>CAL. CHECK<br>AIR BLANK |  | BRAC VALUE<br>.000<br>.000<br>.000 |  | TIME<br>13:44<br>13:45<br>13:45       |  |
| TIME FIRST OBSERVED  |  | SUBJECT'S NAME                               |  | LAB                                |  | INSTRUMENT LOCATION                   |  |
| 461160   |  | OPERATOR                                     |  | 461160                             |  | ADDITIONAL INFORMATION AND/OR REMARKS |  |
| INTOXILYZER® INSTRUMENT PRINTER CARD   |  | © 1986 by CMI INC.                           |  | CMI                                |  |                                       |  |

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT P

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 5000\*

STANDARD CALIBRATION CHECK PROCEDURE *sh*

Agency Tucson Police Department Date 12/2/05 Time 1330  
Intoxilyzer Serial # 166-005004 Location LAB  
QA Specialist Jennifer Coombs 46160  
(Print Name)

- (✓) 1. Pour a standard alcohol solution of known value into a clean dry simulator jar and assemble the simulator. Insure that a tight seal has been made. Standard value: 0.080 AC
- (✓) 2. Turn on the simulator and allow the temperature to reach 34°C ± .2°C.
- (✓) 3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- (✓) 4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- ✓✓ 5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- ✓✓ 6. Push Start Test button or press enter on keyboard.
- ✓✓ 7. Insert card in response to display.
- ✓✓ 8. Air blank completed.
- ✓✓ 9. Calibration check completed. Test results 0.080 AC & 0.079 AC
- ✓✓ 10. Air blank completed.
- ✓✓ 11. When display reads Test Complete remove evidence card. Attach the card to the completed checklist.
- (✓) 12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE

Jennifer L. Coombs

\*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DHS/DSLS/Form C116 (Rev. 7-93)

TUCSON POLICE DEPARTMENT  
INTOXILYZER - ALCOHOL ANALYZER  
AZ MODEL 5000 SN 66-005224  
12/02/05 13:50

TEST BRAC VALUE TIME  
AIR BLANK .000 13:49  
INVALID SAMPLE .004 13:50  
AIR BLANK .000 13:50

BEGIN NEW DEPRIVATION PERIOD

SUBJECT'S NAME

TIME FIRST OBSERVED

LAB

INSTRUMENT LOCATION

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

Sweet Breath

INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

CMI

TUCSON POLICE DEPARTMENT  
INTOXILYZER - ALCOHOL ANALYZER  
AZ MODEL 5000 SN 66-005224  
12/02/05 13:54

TEST BRAC VALUE TIME  
AIR BLANK .000 13:51  
\*SUBJECT TEST .075 13:54  
AIR BLANK .000 13:54

\* DEFICIENT SAMPLE - VALUE  
PRINTED WAS HIGHEST OBTAINED.

SUBJECT'S NAME

TIME FIRST OBSERVED

LAB

INSTRUMENT LOCATION

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

CMI

6 of 6  
file



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

INTOXILYZER SERVICE RECORD *Bh*

*1 of 1*  
*file*

Date: December 2, 2005

QAS: Jennifer Coombs 466160

Intoxilyzer Serial Number: 66-005004

Reason For Service: The power switch was broken by the officer.

Date Out of Service: 12/2/05

Calibration Check completed prior to removal from service

Date: 12/2/05

SQAP completed prior to removal from service

Date: 12/2/05

Corrective Action:

☒ Sent to Applied Electronics