



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000 332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

T54
5-2-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4-28-11 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 3-31-11 and 4-28-11
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


4-28-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 28, day of April, 20 11.

Isabel C. Trujillo

Notary Public

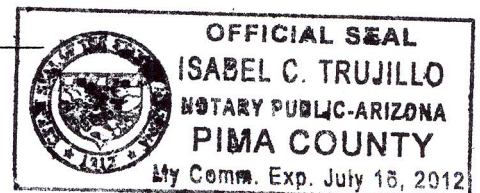

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1007
cm

TS4
5.2.2011

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 4-28-11 TIME 1154

INTOXILYZER SERIAL # 80-000 332 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS None.

SIGNATURE Clare Conley

DPS Form Exh G-4 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
04/28/2011 11:56:35

2067
CH

Standard Lot#: 719919
Last Changed By: W. HONOMICHL #42992

TS4

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:57:22
Diagnostic Test	Pass	11:57:50
Air Blank	0.000	11:58:17
0.100 Cal Check	0.098	11:58:36
Air Blank	0.000	11:59:05
Subject Test	0.000	12:00:11
Air Blank	ABT*	12:00:18
Air Blank	0.000	12:00:46

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
04/28/2011 12:00:50

3067
CM
134

Standard Lot#: 719919
Last Changed By: W. HONOMICHL #42992

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	12:01:13
Air Blank	IPS**	12:01:25

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000332

04/28/2011

12:01:30

Standard Lot#: 719919

Last Changed By: W. HONOMICHL #42992

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:02:16
Diagnostic Test	Pass	12:02:43
Air Blank	0.000	12:03:11
0.100 Cal Check	0.097	12:03:29
Air Blank	0.000	12:03:59
Subject Test	0.078*	12:07:07
Air Blank	0.000	12:07:40
Air Blank	0.000	12:08:07
Subject Test	XXX**	12:08:24
Air Blank	0.000	12:08:58

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

04/28/2011

12:09:16


Standard Lot#: 719919

Last Changed By: W. HONOMICHL #42992

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:09:44
Cal Check	0.097	12:10:02
Air Blank	0.000	12:10:32

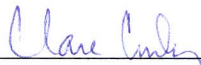

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

677
CH

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSY
5.2.2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 4-28-11 TIME 1154
Intoxilyzer Serial # 80-000 332 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

04/28/2011

11:54:52

7-67
CN

Standard Lot#: 719919

Last Changed By: W. HONOMICHL #42992

154

QAS: C. CONLEY #34517

TUCSON PD

Test g/210L Time

Air Blank 0.000 11:55:19
Cal Check 0.097 11:55:38
Air Blank 0.000 11:56:08

Clare Conley

Operator's Signature