



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

734

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12/29/06 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 12/1/06 and 12/29/06 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Jennifer L. Coombs  
Jennifer L. Coombs

1/4/07

ADPS QAS Permit # 33120

SUBSCRIBED AND SWORN to before me this 4, day of January, 2007.

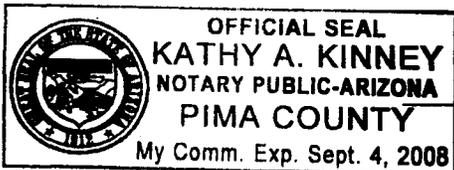

Kathy A. Kinney  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

1 of 8  
JL

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE **74**

QA SPECIALIST Jennifer Coombs 96160 AGENCY Tucson Police Department

DATE December 29, 2006 TIME 0911

INTOXILYZER SERIAL # 80-000323 LOCATION ODD

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed ✓✓✓
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed ✓✓  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Jennifer L. Coombs

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

12/29/2006

09:42:31

*PH*

*2 of 8  
file*

-----  
Standard Lot#: 502602I

Last Changed By: J. COOMBS #46160

-----  
QAS: J. COOMBS #46160

TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	09:42:59
Cal Check	0.097	09:43:17
Air Blank	0.000	09:43:47

*Jennifer L. Coombs*  
-----  
Operator's Signature

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000323  
12/29/2006

09:20:50

*BJ*

*3 of 8  
file*

---

Standard Lot#: 502602I  
Last Changed By: J. COOMBS #46160

---

QAS: J. COOMBS #46160  
TUCSON PD CRIME LAB

---

Operator: J. COOMBS #46160  
TUCSON PD CRIME LAB

---

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:22:10
Diagnostic Test	Pass	09:22:37
Air Blank	0.000	09:23:05
0.100 Cal Check	0.098	09:23:23
Air Blank	0.000	09:23:52
Subject Test	0.000	09:24:14
Air Blank	0.000	09:24:42
Five Minute Wait	ABT*	
Air Blank	0.000	09:26:39

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000323  
12/29/2006

09:28:59

W

4 of 8  
file

---

Standard Lot#: 502602I  
Last Changed By: J. COOMBS #46160

---

QAS: J. COOMBS #46160  
TUCSON PD CRIME LAB

---

Operator: J. COOMBS #46160  
TUCSON PD CRIME LAB

---

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	IPS*	09:29:52
Air Blank	0.000	09:30:21

\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000323  
12/29/2006

09:30:25

76

5 of 8  
me

---

Standard Lot#: 502602I  
Last Changed By: J. COOMBS #46160

---

QAS: J. COOMBS #46160  
TUCSON PD CRIME LAB

---

Operator: J. COOMBS #46160  
TUCSON PD CRIME LAB

---

Subject: SQAP, TEST

DOB: 11/11/1111

Sex: F

Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:31:26
Diagnostic Test	Pass	09:31:53
Air Blank	0.000	09:32:21
0.100 Cal Check	RFI*	09:32:36
Air Blank	0.000	09:33:07

\*RFI Detect

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000323  
12/29/2006

09:33:25

TS4

leaf 8  
file

---

Standard Lot#: 502602I  
Last Changed By: J. COOMBS #46160

---

QAS: J. COOMBS #46160  
TUCSON PD CRIME LAB

---

Operator: J. COOMBS #46160  
TUCSON PD CRIME LAB

---

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:34:18
Diagnostic Test	Pass	09:34:46
Air Blank	0.000	09:35:13
0.100 Cal Check	0.097	09:35:32
Air Blank	0.000	09:36:01
Subject Test	0.040*	09:39:10
Air Blank	0.000	09:39:41
Air Blank	0.000	09:40:08
Subject Test	XXX**	09:41:00
Air Blank	0.000	09:41:33

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

EXHIBIT G-2

7 of 8  
file

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

PH

QA Specialist Jennifer Coombs 46160 AGENCY Tucson Police Department  
DATE December 29, 2006 TIME 0911  
Intoxilyzer Serial # 80-000323 Location 000

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C + 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Jennifer L. Coombs

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
12/29/2006 09:18:05


8 of 8  
JL

P4

-----  
Standard Lot#: 502602I  
Last Changed By: J. COOMBS #46160  
-----

QAS: J. COOMBS #46160  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	09:18:32
Cal Check	0.097	09:18:51
Air Blank	0.000	09:19:20

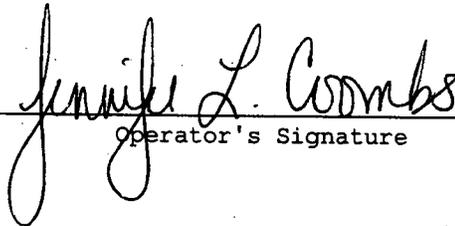  
\_\_\_\_\_  
Operator's Signature