



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

me

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2.22.2007 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 1.29.2007 and 2.22.2007 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on                     .  
The instrument was returned to service on                     .

Terry S. Gallegos  
Terry S. Gallegos

March 12, 2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 12, day of March, 2007.

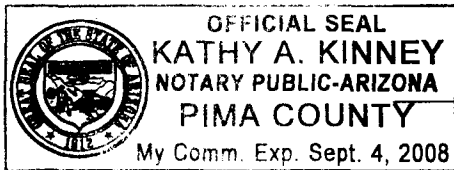

Kathy A Kinney  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1089  
TSY mm

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 2.22.2007 TIME 1649

INTOXILYZER SERIAL # 80-000323 LOCATION Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed ✓✓✓  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ✓ NO       

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Terry Gallegos

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000323  
02/22/2007

16:51:21

2059  
TSh<sup>no</sup>

Standard Lot#: 627101G  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
~~AZ-DPS~~ <sup>TSh</sup> POLICE DEPARTMENT/CRIME LAB

Operator: T. GALLEGOS #28847  
~~AZ-DPS~~ <sup>TSh</sup> POLICE DEPARTMENT/CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:52:18
Diagnostic Test	Pass	16:52:45
Air Blank	0.000	16:53:13
0.100 Cal Check	0.100	16:53:32
Air Blank	0.000	16:54:01
Subject Test	0.000	16:54:21
Air Blank	0.000	16:54:48
Five Minute Wait	ABT*	
Air Blank	IPS**	16:55:03

\*Sequence Aborted

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000323  
02/22/2007

16:58:59

3059  
TSU m

Standard Lot#: 627101G  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
~~AZ DPS~~ TSU POLICE DEPARTMENT/CRIME LAB

Operator: T. GALLEGOS #28847  
~~AZ DPS~~ TSU POLICE DEPARTMENT/CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:59:54
Diagnostic Test	Pass	17:00:21
Air Blank	0.000	17:00:48
0.100 Cal Check	0.099	17:01:07
Air Blank	0.000	17:01:36
Subject Test	RFI*	17:01:57
Air Blank	0.000	17:02:26

\*RFI Detect

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000323  
02/22/2007

17:05:45

Standard Lot#: 627101G  
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847  
~~AZ DES~~ **POLICE DEPARTMENT/CRIME LAB**

Operator: ~~BY~~ T. GALLEGOS #28847  
TPD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:07:02
Diagnostic Test	Pass	17:07:29
Air Blank	0.000	17:07:56
0.100 Cal Check	0.100	17:08:15
Air Blank	0.000	17:08:44
Subject Test	0.086*	17:11:53
Air Blank	0.000	17:12:25
Air Blank	0.000	17:12:52
Subject Test	XXX**	17:13:17
Air Blank	0.000	17:13:52

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

02/22/2007

17:14:15

Standard Lot#: 627101G

Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847

~~AZ DPS~~

**POLICE DEPARTMENT/CRIME LAB**

Test	g/210L	Time
Air Blank	0.000	17:14:42
Cal Check	0.099	17:15:01
Air Blank	0.000	17:15:31

  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 2.22.2007 TIME 1649  
Intoxilyzer Serial # 80-000 323 Location Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

02/22/2007

16:49:27

705 9  
BY ms

-----  
Standard Lot#: 627101G

Last Changed By: T. GALLEGOS #28847

-----  
QAS: T. GALLEGOS #28847

~~AZ DPS~~

**TS4 POLICE DEPARTMENT/CRIME LAB**

-----  
Test g/210L Time  
-----  
Air Blank 0.000 16:49:54  
Cal Check 0.099 16:50:13  
Air Blank 0.000 16:50:42

  
Operator's Signature





TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

8059  
B4 m

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

Date: 2.22.2007  
QAS: Terry Gallegos 28947  
Intoxilyzer Serial Number: 80-000323  
New Lot Number: 6271019  
New Expiration Date: 10.03.2008

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
Enter Target Value: 0.100  
Enter Control Lot #: From the tank label  
Enter Exp Date: From the tank label (mm/dd/yyyy)  
Edit info Y/N? Y or N  
Update Std Info? Y ENTER  
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

02/22/2007

16:43:50

Standard Lot#: 627101G

Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847

~~AZ DPS~~

TS4

POLICE DEPARTMENT/CRIME LAB

Test	g/210L	Time
Air Blank	0.000	16:44:18
Cal Check	0.100	16:44:37
Air Blank	0.000	16:45:06

  
Operator's Signature

908 9  
TS4 mo