



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

B4
7-10-2007

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7-5-07 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 6-7-07 and 7-5-07 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

7-9-07

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 9, day of July, 2007.

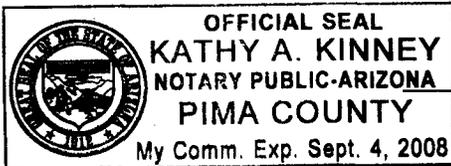

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

β4

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department
DATE 7-5-07 TIME 1409
INTOXILYZER SERIAL # 80-000323 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 10 cu

B4

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000323

07/05/2007

14:13:14

Standard Lot#: 706602I

Last Changed By: C. CONLEY

#34517

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:14:14
Diagnostic Test	Pass	14:14:42
Air Blank	0.000	14:15:09
0.100 Cal Check	0.099	14:15:28
Air Blank	0.000	14:15:57
Subject Test	0.000	14:16:18
Air Blank	0.000	14:16:45
Five Minute Wait	ABT*	
Air Blank	0.000	14:17:25

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 10 cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
07/05/2007

14:18:03

Standard Lot#: 706602I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:19:00
Diagnostic Test	Pass	14:19:28
Air Blank	0.000	14:19:55
0.100 Cal Check	RFI*	14:20:07
Air Blank	0.000	14:20:38

*RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 10 Cn

74

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
07/05/2007 14:21:12

Standard Lot#: 706602I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	14:21:47
Air Blank	0.000	14:22:16

*Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5 of 10 cu

BA

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
07/05/2007 14:23:21

Standard Lot#: 706602I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:24:20
Diagnostic Test	Pass	14:24:48
Air Blank	0.000	14:25:15
0.100 Cal Check	0.098	14:25:34
Air Blank	0.000	14:26:03
Subject Test	0.096*	14:29:11
Air Blank	0.000	14:29:44
Air Blank	0.000	14:30:12
Subject Test	XXX**	14:30:35
Air Blank	0.000	14:31:09

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

6 of 10 cu

Bh

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
07/05/2007 14:33:49

Standard Lot#: 706602I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:34:17
Cal Check	0.098	14:34:35
Air Blank	0.000	14:35:05

Clare Conley
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

7 of 10 cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

Bh

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 7-5-07 TIME 1409
Intoxilyzer Serial # 80-000323 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)


8 of 10 cu
B4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
07/05/2007 14:08:44

Standard Lot#: 706602I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:09:12
Cal Check	0.099	14:09:30
Air Blank	0.000	14:10:00

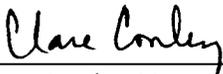

Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

9 of 10 cu

134

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

Date: 7-5-07
QAS: Clare Conley
Intoxilyzer Serial Number: 80-000323
New Lot Number: 706602I
New Expiration Date: 11 Mar 2009

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change." Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing." Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0.100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE

Clare Conley

10 of 10 cu


734

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
07/05/2007 14:07:09

Standard Lot#: 706602I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:07:36
Cal Check	0.100	14:07:55
Air Blank	0.000	14:08:24

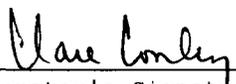

Operator's Signature