



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

TJH  
5-2-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5-2-08 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 4-4-08 and 5-2-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

5-2-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 2, day of May, 2008.

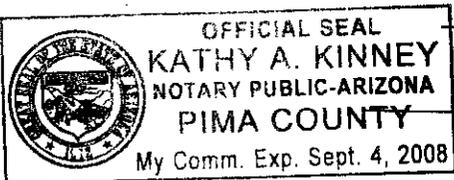

Kathy A. Kinney  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

568  
cu  
BY  
5.2-08

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 5-2-08 TIME 1134

INTOXILYZER SERIAL # 80-000323 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000323  
05/02/2008

11:35:56

2 of 8 cu  
TS4

Standard Lot#: 706601I  
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:36:43
Diagnostic Test	Pass	11:37:10
Air Blank	0.000	11:37:37
0.100 Cal Check	0.096	11:37:56
Air Blank	0.000	11:38:26
Subject Test	0.000	11:39:02
Air Blank	IPS*	11:39:13
Air Blank	0.000	11:39:41

\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000323  
05/02/2008

11:39:46

3.8 cu  
BY

Standard Lot#: 706601I  
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	11:40:17
Air Blank	0.000	11:40:46

\*RFI Detect

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000323  
05/02/2008

11:41:03

4/28/08  
TS4

Standard Lot#: 706601I  
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	ABT*	11:41:38
Air Blank	0.000	11:42:06

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000323  
05/02/2008

11:42:09

5 of 8 cu  
TS4

Standard Lot#: 706601I  
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:42:59
Diagnostic Test	Pass	11:43:26
Air Blank	0.000	11:43:53
0.100 Cal Check	0.097	11:44:12
Air Blank	0.000	11:44:42
Subject Test	0.237*	11:47:50
Air Blank	0.000	11:48:23
Air Blank	0.000	11:48:50
Subject Test	XXX**	11:49:18
Air Blank	0.000	11:49:52

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

05/02/2008

11:56:45

6 of 8 cm

B4

-----  
Standard Lot#: 706601I

Last Changed By: S. CIONE #50616  
-----

QAS: C. CONLEY #34517

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:57:12
Cal Check	0.096	11:57:31
Air Blank	0.000	11:58:00

*Clare Conley*

\_\_\_\_\_  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

7 of 8 cu

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

B4  
5-2-08

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 5-2-08 TIME 1134  
Intoxilyzer Serial # 80-000323 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
05/02/2008

11:34:26

-----  
Standard Lot#: 706601I  
Last Changed By: S. CIONE #50616  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:34:53
Cal Check	0.097	11:35:12
Air Blank	0.000	11:35:42

*Clare Conley*

Operator's Signature

*8 of 8 cu*

*TS4*