



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

134
2-17-09

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/6/09 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 1/9/09 and 2/6/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


2/17/09

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 17th, day of February, 2009.

A.M. Hernández

Notary Public

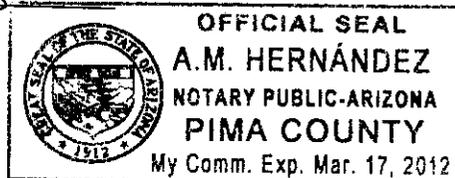

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

B4
2.17.09

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 2/6/09 TIME 13:31

INTOXILYZER SERIAL # 80-000323 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7sc

134

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
02/06/2009

13:33:18

Standard Lot#: 719802I
Last Changed By: J. TRAINING #99999

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:34:49
Diagnostic Test	Pass	13:35:16
Air Blank	0.000	13:35:44
0.100 Cal Check	0.099	13:36:03
Air Blank	0.000	13:36:32
Subject Test	0.000	13:37:01
Air Blank	0.000	13:37:29
Five Minute Wait	ABT*	
Air Blank	0.000	13:38:06

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 7 sc

134

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
02/06/2009

13:38:33

Standard Lot#: 719802I
Last Changed By: J. TRAINING #99999

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	13:39:03
Air Blank	IPS**	13:39:11

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
02/06/2009

13:40:42

4 of 7 sc

134

Standard Lot#: 719802I
Last Changed By: J. TRAINING #99999

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:41:45
Diagnostic Test	Pass	13:42:12
Air Blank	0.000	13:42:40
0.100 Cal Check	0.099	13:42:58
Air Blank	0.000	13:43:28
Subject Test	0.064*	13:46:36
Air Blank	0.000	13:47:08
Air Blank	0.000	13:47:35
Subject Test	XXX**	13:48:02
Air Blank	0.000	13:48:35

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 7 SC

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
02/06/2009 13:49:10

Standard Lot#: 719802I
Last Changed By: J. TRAINING #99999

BY

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:49:38
Cal Check	0.097	13:49:57
Air Blank	0.000	13:50:26

Sarah Cione
Operator's Signature

EXHIBIT G-2

6 of 7 SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
2-17-09

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 2/6/09 TIME 13:31
Intoxilyzer Serial # 80-000323 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- () 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
02/06/2009 13:31:32

7 of 7sc

154

Standard Lot#: 719802I
Last Changed By: J. TRAINING #99999

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:32:00
Cal Check	0.098	13:32:19
Air Blank	0.000	13:32:48

Sarah Cione
Operator's Signature