



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000 323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

By 5.11.09

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5-1-09 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 4-2-09 and 5-1-09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A. The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


May 8, 2009

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 8th day of May, 2009.

A. Hernandez

Notary Public

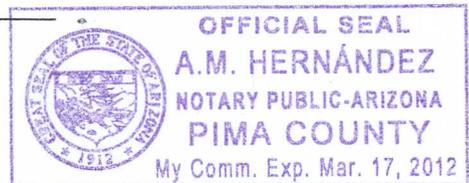

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1087
eh

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

By
5.11.09

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 5-1-09 TIME 1435

INTOXILYZER SERIAL # 80-000323 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2087
CU
184

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
05/01/2009 14:37:22

Standard Lot#: 811201I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:38:10
Diagnostic Test	Pass	14:38:38
Air Blank	0.000	14:39:05
0.100 Cal Check	0.099	14:39:24
Air Blank	0.000	14:39:53
Subject Test	0.000	14:40:10
Air Blank	ABT*	14:40:16
Air Blank	0.000	14:40:45

*Sequence Aborted

Not a Successfully Completed Test Sequence

3587
Ch
184

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
05/01/2009 14:40:48

Standard Lot#: 811201I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	14:41:22
Air Blank	IPS**	14:41:33

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4567
ch
PS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
05/01/2009 14:41:40

Standard Lot#: 811201I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:42:31
Diagnostic Test	Pass	14:42:59
Air Blank	0.000	14:43:26
0.100 Cal Check	0.098	14:43:45
Air Blank	0.000	14:44:14
Subject Test	0.214*	14:47:23
Air Blank	0.000	14:47:56
Air Blank	0.000	14:48:23
Subject Test	XXX**	14:48:52
Air Blank	0.000	14:49:25

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
05/01/2009 14:49:39

5567
lu
B4

Standard Lot#: 811201I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:50:06
Cal Check	0.097	14:50:25
Air Blank	0.000	14:50:55

Clare Conley

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6587
Cm
Joh 5.11.09

QA Specialist Clara Conley AGENCY Tucson Police Department
DATE 5-1-09 TIME 1435
Intoxilyzer Serial # 80-000323 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clara Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
05/01/2009 14:35:33

7567

CH

Standard Lot#: 811201I
Last Changed By: S. CIONE #50616

PH

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:36:00
Cal Check	0.098	14:36:19
Air Blank	0.000	14:36:48

Clare Conley

Operator's Signature