



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/29/09 Cal 1 0.096 AC Cal 2 0.096 AC

Between the dates of 5/28/09 and 6/29/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


6/29/09

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 29th day of June, 20 09.

A.M. Hernandez

Notary Public

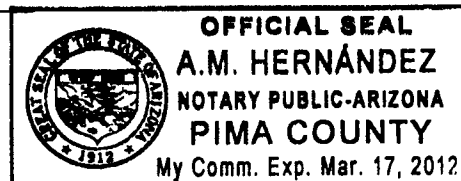

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

T54
6-29-09

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 6/29/09 TIME 10:20
INTOXILYZER SERIAL # 80-000323 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 096 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE

Sarah Cione

2 of 7sc

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000323

06/29/2009

10:22:12

104

Standard Lot#: 811201I

Last Changed By: C. CONLEY

#34517

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:23:01
Diagnostic Test	Pass	10:23:29
Air Blank	0.000	10:23:56
0.100 Cal Check	0.096	10:24:15
Air Blank	0.000	10:24:44
Subject Test	0.000	10:25:00
Air Blank	0.000	10:25:27
Five Minute Wait	ABT*	
Air Blank	0.000	10:26:08

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7sc

164

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000323

06/29/2009

10:26:19

Standard Lot#: 811201I

Last Changed By: C. CONLEY

#34517

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP,TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:26:46
Air Blank	IPS**	10:26:54

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000323

06/29/2009

10:27:28

154

Standard Lot#: 811201I

Last Changed By: C. CONLEY

#34517

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:28:36
Diagnostic Test	Pass	10:29:03
Air Blank	0.000	10:29:31
0.100 Cal Check	0.097	10:29:49
Air Blank	0.000	10:30:19
Subject Test	0.058*	10:33:27
Air Blank	0.000	10:33:59
Air Blank	0.000	10:34:27
Subject Test	XXX**	10:34:52
Air Blank	0.000	10:35:25

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 7 sc

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

06/29/2009

10:36:12

104

Standard Lot#: 811201I

Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:36:40
Cal Check	0.096	10:36:58
Air Blank	0.000	10:37:28

Sarah Cione

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSY
6-29-09

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 6/29/09 TIME 10:20
Intoxilyzer Serial # 80-000323 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

06/29/2009

10:20:32


Standard Lot#: 811201I

Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:20:59
Cal Check	0.096	10:21:18
Air Blank	0.000	10:21:48

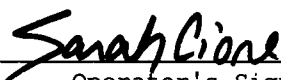

Operator's Signature