



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

154
2-16-2010

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 02/11/10 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 01/13/10 and 02/11/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
 The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub

02/11/10

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 11th day of February, 2010.

A.M. Hernandez

Notary Public



*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

197
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

TS
2-16-2010

QA SPECIALIST CHERYL STRAUB AGENCY Tucson Police Department
DATE 02/11/10 TIME 11:39
INTOXILYZER SERIAL # 80-000323 LOCATION CRIME LAB

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Cheryl Straub

217
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
02/11/2010 11:41:07

124

Standard Lot#: 811201I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:42:01
Diagnostic Test	Pass	11:42:29
Air Blank	0.000	11:42:56
0.100 Cal Check	0.099	11:43:15
Air Blank	0.000	11:43:44
Subject Test	0.000	11:44:03
Air Blank	ABT*	11:44:15
Air Blank	0.000	11:44:43

*Sequence Aborted

Not a Successfully Completed Test Sequence

377
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
02/11/2010 11:45:09

TS

Standard Lot#: 811201I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	11:45:55
Air Blank	IPS**	11:46:07

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

417
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
02/11/2010 11:46:44

154

Standard Lot#: 8112011
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB* #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:47:37
Diagnostic Test	Pass	11:48:05
Air Blank	0.000	11:48:32
0.100 Cal Check	0.098	11:48:51
Air Blank	0.000	11:49:20
Subject Test	0.106*	11:52:29
Air Blank	0.000	11:53:01
Air Blank	0.000	11:53:28
Subject Test	XXX**	11:54:02
Air Blank	0.000	11:54:35

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

577
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
02/11/2010 11:55:13

Standard Lot#: 8112011
Last Changed By: C. CONLEY #34517

104

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:55:40
Cal Check	0.097	11:55:59
Air Blank	0.000	11:56:28



Operator's Signature

EXHIBIT G-2

677
CSTHIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

154
2-16-2010

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
 DATE 02/11/10 TIME 11:39
 Intoxilyzer Serial # 80-000323 Location CRIME/LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
 OR
 Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

777
CS

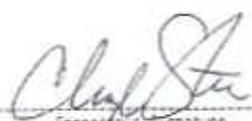
TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-100323
02/11/2010 11:39:09

Standard Lot#: 8112011
Last Changed By: C. CONLEY #34517

DRS: C. STRUB #58772
TUCSON PD

Test	g/210L	Time
Rt Blank	0.000	11:39:36
Cal Check	0.098	11:39:55
Rt Blank	0.000	11:40:24


Operator's Signature