



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

By
3.11.2010

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 03/11/10 Cal 1 0.099 AC Cal 2 0.097 AC

Between the dates of 02/11/10 and 03/11/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


03/11/10

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 11th, day of March, 2010.

A.M. Hernandez

Notary Public

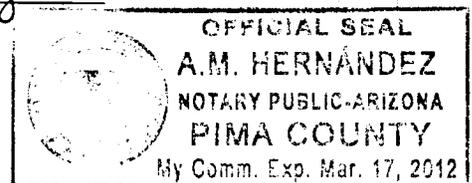

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 7
CS

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

TS4
3-11-2010

QA SPECIALIST Cheryl Straub AGENCY Tucson Police Department
DATE 03/11/10 TIME 10:09
INTOXILYZER SERIAL # 80-000323 LOCATION CRIME LAB

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

277
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
03/11/2010

10:11:17

TS4

Standard Lot#: 811201I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:12:09
Diagnostic Test	Pass	10:12:36
Air Blank	0.000	10:13:04
0.100 Cal Check	0.098	10:13:22
Air Blank	0.000	10:13:52
Subject Test	0.000	10:14:11
Air Blank	ABT*	10:14:18
Air Blank	0.000	10:14:46

*Sequence Aborted

Not a Successfully Completed Test Sequence

377
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
03/11/2010 10:15:03

TS4

Standard Lot#: 811201I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:15:35
Air Blank	IPS**	10:15:43

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

477
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
03/11/2010

10:15:49

DS

Standard Lot#: 811201I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:16:36
Diagnostic Test	Pass	10:17:04
Air Blank	0.000	10:17:31
0.100 Cal Check	0.099	10:17:50
Air Blank	0.000	10:18:19
Subject Test	0.226*	10:21:27
Air Blank	0.000	10:22:00
Air Blank	0.000	10:22:27
Subject Test	XXX**	10:22:48
Air Blank	0.000	10:23:22

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

577
C8

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

03/11/2010

10:24:32

Standard Lot#: 811201I

Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772

TUCSON PD

TS4

Test	g/210L	Time
Air Blank	0.000	10:24:59
Cal Check	0.097	10:25:18
Air Blank	0.000	10:25:47



Operator's Signature

697
es

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TSK
3.11.2010

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 03/11/10 TIME 10:09
Intoxilyzer Serial # 80-000323 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution n/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

777
cs

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

03/11/2010

10:09:40

Standard Lot#: 811201I

Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772

TUCSON PD

154

Test	g/210L	Time
Air Blank	0.000	10:10:07
Cal Check	0.099	10:10:26
Air Blank	0.000	10:10:55



Operator's Signature