



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number

80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

By  
5.11.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 05/05/11 Cal 1 0.094 AC Cal 2 0.094 AC

Between the dates of 04/08/2011 and 05/05/2011 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


05/11/2011

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 11, day of May, 2011.

Ant. Hernandez

Notary Public

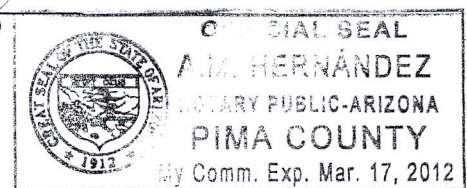

\*Includes documents specified in A.A.C. R13-10-104.A.5

187cs

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

B4  
5.11.2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department  
DATE 05/05/2011 TIME 10:24  
INTOXILYZER SERIAL # 80-000323 LOCATION CRIME LAB

(✓) 1. Display reads "PUSH BUTTON TO START"

## DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
(✓) 2. Date check.

## OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 094 AC

Instrument is operating properly and accurately. YES ✓ NO     

## COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Cheryl Str

287  
CS

154

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
05/05/2011 10:26:28

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SOAP, SOAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:27:24
Diagnostic Test	Pass	10:27:53
Air Blank	0.000	10:28:20
0.100 Cal Check	0.095	10:28:39
Air Blank	0.000	10:29:08
Subject Test	0.000	10:29:54
Air Blank	ABT*	10:30:00
Air Blank	0.000	10:30:28

\*Sequence Aborted

Not a Successfully  
Completed  
Test Sequence

387  
CS

TS4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 60-000323  
05/05/2011 10:31:13

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

OAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SOAP, SOAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:31:46
Air Blank	IPS**	10:31:54

\*RFI Detect

\*\*Improper Sample

Not a Successfully  
Completed  
Test Sequence

487  
CS

B4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
15/05/2011 10:32:32

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

OAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:33:22
Diagnostic Test	Pass	10:33:51
Air Blank	0.000	10:34:19
0.100 Cal Check	0.094	10:34:37
Air Blank	0.000	10:35:07
Subject Test	0.237*	10:36:15
Air Blank	0.000	10:36:47
Air Blank	0.000	10:39:15
Subject Test	XXX**	10:39:34
Air Blank	0.000	10:40:08

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

Not a Successfully  
Completed  
Test Sequence

507  
CS

PM

INTOXILYZER 6000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-100323  
05/05/2011 10:40:35

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

OAS: C. STRAUB #50772  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:41:03
Cal Check	0.094	10:41:21
Air Blank	0.000	10:41:51

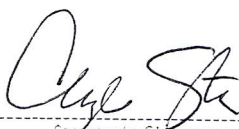
  
Operator's Signature



EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department  
DATE 05/05/2011 TIME 10:24  
Intoxilyzer Serial # 80-000323 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 094 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

787  
CS

TS4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
05/05/2011 10:24:45

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

OAS: C. STRAUB #50772  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:25:13
Cal Check	0.094	10:25:32
Air Blank	0.000	10:26:01

  
Operator's Signature