



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
7-5-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 06/30/11 Cal 1 0.093 AC Cal 2 0.091 AC

Between the dates of 06/02/11 and 06/30/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

*Cheryl K. Straub*


Cheryl K. Straub

06/30/11

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 30th day of June, 2011.

*Isabel C. Trujillo*  
Notary Public

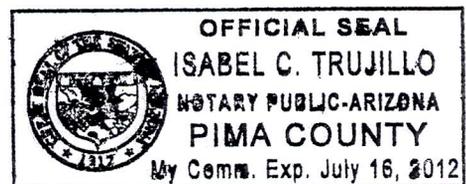

\*Includes documents specified in A.A.C. R13-10-104.A.5

187  
CS

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES

INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TB4  
7-5-2011

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department

DATE 04/30/2011 TIME 08:46

INTOXILYZER SERIAL # 80-000323 LOCATION CRIME LAB

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

( ✓ ) 1. Clock time check.

( ✓ ) 2. Date check.

OPERATIONAL TESTS

( ✓ ) 1. Alcohol-free subject test result 0. 060 AC.

( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed

( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.

( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 091 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Cheryl Straub

287  
CS

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000323  
06/30/2011

08:47:58

TS4

Standard Lot#: 642845  
Last Changed By: T. GALLEGOS #28847

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:48:48
Diagnostic Test	Pass	08:49:15
Air Blank	0.000	08:49:43
0.100 Cal Check	0.092	08:50:01
Air Blank	0.000	08:50:31
Subject Test	0.000	08:51:03
Air Blank	ABT*	08:51:10
Air Blank	0.000	08:51:38

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

387  
CS

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000323  
06/30/2011 08:51:42

TS4

Standard Lot#: 642845  
Last Changed By: T. GALLEGOS #28847

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:52:11
Air Blank	IPS**	08:52:19

\*RFI Detect  
\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

477  
CS

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000323  
06/30/2011 08:52:26

TS4

Standard Lot#: 642845  
Last Changed By: T. GALLEGOS #28847

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:53:30
Diagnostic Test	Pass	08:53:58
Air Blank	0.000	08:54:25
0.100 Cal Check	0.093	08:54:44
Air Blank	0.000	08:55:13
Subject Test	0.140*	08:58:22
Air Blank	0.000	08:58:54
Air Blank	0.000	08:59:21
Subject Test	XXX**	08:59:43
Air Blank	0.000	09:00:17

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

587  
CS

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
06/30/2011 09:00:33

TS4

-----  
Standard Lot#: 642845  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:01:01
Cal Check	0.091	09:01:19
Air Blank	0.000	09:01:49



Operator's Signature

637  
CS

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

TB4  
7.5.2011

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department  
DATE 06/30/2011 TIME 08:46  
Intoxilyzer Serial # 80-000323 Location CRIME LAB

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 093 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

777  
CS


INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
06/30/2011 08:46:23

TS4

-----  
Standard Lot#: 642845  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:46:50
Cal Check	0.093	08:47:09
Air Blank	0.000	08:47:38

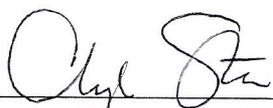  
-----  
Operator's Signature