



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
12.22.2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80 - 000 323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-15-11 Cal 1 0.103 AC Cal 2 0.101 AC

Between the dates of 11-18-11 and 12-15-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


12-22-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 22nd, day of December, 20 11.

Diane Quiroz

Notary Public

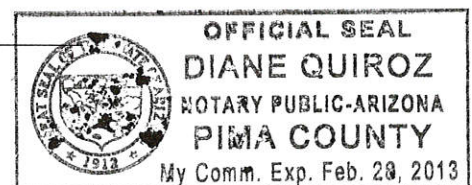

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department
DATE 12-15-11 TIME 1355
INTOXILYZER SERIAL # 80-000323 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

None.

SIGNATURE

Clare Conley

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
12/15/2011

13:57:19

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1911
Sex: M Weight: 130

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:58:36
Diagnostic Test	Pass	13:59:10
Air Blank	0.000	13:59:37
0.100 Cal Check	0.102	13:59:58
Air Blank	0.000	14:00:28
Subject Test	0.000	14:00:44
Air Blank	ABT*	14:00:51
Air Blank	RFI**	14:00:58

*Sequence Aborted

**RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

Core Version: 8105.48

12/15/2011

14:01:08

Standard Lot#: 20910100A1

Last Changed By: S. CIONE

#50616

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1911

Sex: M

Weight: 130

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	14:01:56
Air Blank	0.000	14:02:25

*Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

Core Version: 8105.48

12/15/2011

14:02:42

Standard Lot#: 20910100A1

Last Changed By: S. CIONE

#50616

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TPD

Subject: TEST, SQAP

DOB: 11/11/1911

Sex: M

Weight: 130

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:03:50
Diagnostic Test	Pass	14:04:24
Air Blank	0.000	14:04:52
0.100 Cal Check	0.102	14:05:12
Air Blank	0.000	14:05:42
Subject Test	0.372*	14:08:49
Air Blank	0.000	14:09:25
Air Blank	0.000	14:09:53
Subject Test	RGE**	14:10:18
Air Blank	0.000	14:10:54

*Deficient Sample

**Range Exceeded

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 8
CH
134

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
12/15/2011

14:11:11

528
CM
JBY

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TPD

Subject: TEST, SQAP
DOB: 11/11/1911
Sex: M Weight: 130

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:12:16
Diagnostic Test	Pass	14:12:50
Air Blank	0.000	14:13:18
0.100 Cal Check	0.101	14:13:39
Air Blank	0.000	14:14:09
Subject Test	XXX*	14:14:26
Air Blank	RFI**	14:14:39

*Invalid Sample

Begin new deprivation period

**RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

Core Version: 8105.48

12/15/2011

14:15:23

Standard Lot#: 20910100A1

Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517

TUCSON PD

Operator: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:15:54
0.100 Cal Check	0.101	14:16:15
Air Blank	0.000	14:16:45

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 12-15-11 TIME 1355
Intoxilyzer Serial # 80-000323 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 103 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

Core Version: 8105.48

12/15/2011

13:55:48

Standard Lot#: 20910100A1

Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517

TUCSON PD

Operator: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:56:19
0.100 Cal Check	0.103	13:56:40
Air Blank	0.000	13:57:10

Preventative Maintenance
Performed -- Timer Reset