



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

BS4
1-18-2012

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/11/12 Cal 1 0.103 AC Cal 2 0.101 AC

Between the dates of 12/15/11 and 1/11/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


1/11/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 11, day of January, 20 12.

A.M. Hernandez

Notary Public

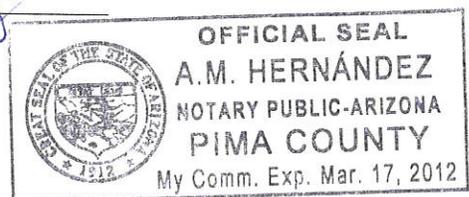

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 7 sc

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

TS4
1-18-2012

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 1/11/12 TIME 1429
INTOXILYZER SERIAL # 80-000323 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7 SC

Sh

INTOXILYZER 8000
 Location: TUCSON POLICE DEPT
 Serial Number: 80-000323
 Core Version: 8105.48
 01/11/2012 14:30:47

Standard Lot#: 20910100A1
 Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
 TUCSON PD

Operator: S. CIONE #50616
 TUCSON PD

Subject: SQAP, SQAP
 DOB: 11/11/1980
 Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:31:30
Diagnostic Test	Pass	14:32:04
Air Blank	0.000	14:32:32
0.100 Cal Check	0.103	14:32:52
Air Blank	0.000	14:33:22
Subject Test	0.000	14:33:42
Air Blank	0.000	14:34:09
Wait	ABT*	
Air Blank	0.000	14:35:12

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 7_{sc}

TS4

INTOXILYZER 8000
 Location: TUCSON POLICE DEPT
 Serial Number: 80-000323
 Core Version: 8105.48
 01/11/2012 14:35:16

Standard Lot#: 20910100A1
 Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
 TUCSON PD

Operator: S. CIONE #50616
 TUCSON PD

Subject: SQAP, SQAP
 DOB: 11/11/1980
 Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	14:35:54
Air Blank	IPS**	14:36:02

*RFI Detect
 **Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 7 sc

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
01/11/2012

14:36:12

Standard Lot#: 20910100A1
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:36:56
Diagnostic Test	Pass	14:37:30
Air Blank	0.000	14:37:57
0.100 Cal Check	0.103	14:38:18
Air Blank	0.000	14:38:48
Subject Test	0.285*	14:41:55
Air Blank	0.000	14:42:28
Air Blank	0.000	14:42:55
Subject Test	XXX**	14:43:22
Air Blank	0.000	14:43:55

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5 of 7_{sc}

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
01/11/2012 14:44:08

Standard Lot#: 20910100A1
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:44:38
0.100 Cal Check	0.101	14:44:59
Air Blank	0.000	14:45:29

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

6 of 7 se

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

TJG
1.18.2012

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 1/11/12 TIME 1429
Intoxilyzer Serial # 80-000323 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0.100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0.103 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
01/11/2012 14:29:14

7 of 7 sc

154

Standard Lot#: 20910100A1
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:29:43
0.100 Cal Check	0.103	14:30:04
Air Blank	0.000	14:30:33

Preventative Maintenance
Performed -- Timer Reset