



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
2.23.2012

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/20/12 Cal 1 0.103 AC Cal 2 0.102 AC

Between the dates of 1/11/12 and 2/20/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


2/21/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 21<sup>st</sup>, day of February, 20 12.

Elma Escarcega

Notary Public

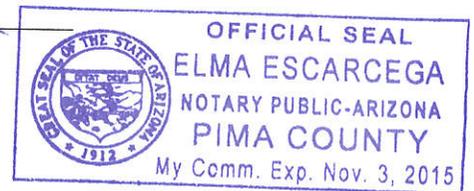

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

1 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

TS4  
2.23.2012

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 2/20/12 TIME 1341

INTOXILYZER SERIAL # 80-000323 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0.100 AC  
Result: 0.102 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS Past due on G-2 by 9 days.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

2 of 7sc  
TJL

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
Core Version: 8105.48  
02/20/2012 13:42:49

---

Standard Lot#: 20910100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:43:40
Diagnostic Test	Pass	13:44:14
Air Blank	0.000	13:44:41
0.100 Cal Check	0.103	13:45:02
Air Blank	0.000	13:45:31
Subject Test	0.000	13:45:50
Air Blank	0.000	13:46:18
Wait	ABT*	
Air Blank	0.000	13:46:52

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

3 of 7 sc

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
Core Version: 8105.48  
02/20/2012 13:47:05

Tsh

---

Standard Lot#: 20910100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	13:47:34
Air Blank	IPS**	13:47:42

\*RFI Detect  
\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
Core Version: 8105.48  
02/20/2012

13:47:46

4 of 7 SC

PS4

---

Standard Lot#: 20910100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:48:29
Diagnostic Test	Pass	13:49:03
Air Blank	0.000	13:49:30
0.100 Cal Check	0.102	13:49:51
Air Blank	0.000	13:50:20
Subject Test	0.035*	13:53:28
Air Blank	0.000	13:54:00
Air Blank	0.000	13:54:27
Subject Test	XXX**	13:54:40
Air Blank	0.000	13:55:14

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
Core Version: 8105.48  
02/20/2012 13:55:28

5 of 7 sc

TS4

-----  
Standard Lot#: 20910100A1  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:55:57
0.100 Cal Check	0.102	13:56:18
Air Blank	0.000	13:56:47

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

6 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TJ  
2-23-2012

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 2/20/12 TIME 1541  
Intoxilyzer Serial # 80-000323 Location Crime Lab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 103 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
Core Version: 8105.48  
02/20/2012 13:41:11

7 of 7sc

Bh

-----  
Standard Lot#: 20910100A1  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:41:40
0.100 Cal Check	0.103	13:42:01
Air Blank	0.000	13:42:30

Preventative Maintenance  
Performed -- Timer Reset