



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4
 4-10-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4/6/12 Cal 1 0.102 AC Cal 2 0.101 AC

Between the dates of 3/12/12 and 4/6/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


4/9/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 9th, day of April, 20 12.

Elma Escarcega

Notary Public

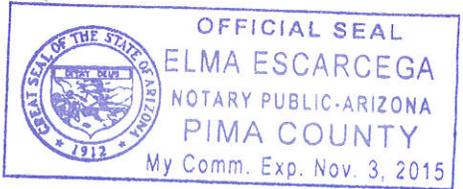

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

1 of 8
SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

T34
4.10.2012

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 4/6/12 TIME 1348

INTOXILYZER SERIAL # 80-000323 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
04/06/2012 13:50:56

2 of 8 SC
BY

Standard Lot#: 20910100A1
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:51:48
Diagnostic Test	Pass	13:52:22
Air Blank	0.000	13:52:49
0.100 Cal Check	0.102	13:53:10
Air Blank	0.000	13:53:40
Subject Test	0.000	13:53:54
Air Blank	0.000	13:54:22
Wait	ABT*	
Air Blank	0.000	13:54:59

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 8 sc

BH

INTOXILYZER 8000
 Location: TUCSON POLICE DEPT
 Serial Number: 80-000323
 Core Version: 8105.48
 04/06/2012 13:55:20

Standard Lot#: 20910100A1
 Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
 TUCSON PD

Operator: S. CIONE #50616
 TUCSON PD

Subject: SQAP, SQAP
 DOB: 11/11/1980
 Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:56:03
Diagnostic Test	Fail*	13:56:07
Air Blank	RFI**	13:56:18

* Voltage/Current Test
 **RFI Detect

→ RFI done during DIAGNOSTICS

Not a Successfully Completed Test Sequence

4 of 8 sc
TSH

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
04/06/2012 14:06:01

Standard Lot#: 20910100A1
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	14:06:25
Air Blank	0.000	14:06:53

*Improper Sample

Not a Successfully Completed Test Sequence

5 of 8_{sc}

TSH

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
04/06/2012 13:56:37

Standard Lot#: 20910100A1
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:57:22
Diagnostic Test	Pass	13:57:56
Air Blank	0.000	13:58:23
0.100 Cal Check	0.102	13:58:44
Air Blank	0.000	13:59:14
Subject Test	0.103*	14:02:21
Air Blank	0.000	14:02:53
Air Blank	0.000	14:03:21
Subject Test	XXX**	14:03:33
Air Blank	0.000	14:04:06

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
04/06/2012 14:04:17

6 of 8_{sc}

B4

Standard Lot#: 20910100A1
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:04:45
0.100 Cal Check	0.101	14:05:06
Air Blank	0.000	14:05:36

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

7 of 8 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

B4
4.10.2012

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department

DATE 4/6/12 TIME 1348

Intoxilyzer Serial # 80-000323 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
04/06/2012 13:48:55

8 of 8 sc

Standard Lot#: 20910100A1
Last Changed By: T. GALLEGOS #28847

BY

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:49:24
0.100 Cal Check	0.102	13:49:45
Air Blank	0.000	13:50:14

Preventative Maintenance
Performed -- Timer Reset