



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

SC
 3/22/13

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3-7-2013 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 2-11-2013 and 3-7-2013 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Terry Gallegos

Terry S. Gallegos


3-13-2013

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 13th day of March, 2013.

Elma Escarcega

Notary Public

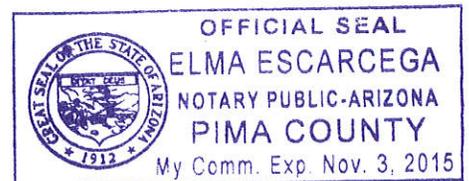

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

LOF7
TBH

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

SC
3/22/13

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos AGENCY Tucson Police Department
DATE 3-7-2013 TIME 1801
INTOXILYZER SERIAL # 80-000323 LOCATION CrimeLab

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Terry Gallegos

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
03/07/2013

18:04:25

2057
TS4

5C

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TPD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:05:45
Diagnostic Test	Pass	18:06:19
Air Blank	0.000	18:06:46
0.100 Cal Check	0.100	18:07:07
Air Blank	0.000	18:07:37
Subject Test	0.000	18:07:55
Air Blank	0.000	18:08:23
Wait	ABT*	
Air Blank	0.000	18:09:00

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
03/07/2013

18:11:05

3 of 7
TS4

SC

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:11:59
Diagnostic Test	Pass	18:12:33
Air Blank	0.000	18:13:01
0.100 Cal Check	RFI*	18:13:15
Air Blank	IPS**	18:13:29

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
03/07/2013

18:13:56

40f7
TSK
SC

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:14:56
Diagnostic Test	Pass	18:15:30
Air Blank	0.000	18:15:57
0.100 Cal Check	0.101	18:16:18
Air Blank	0.000	18:16:48
Subject Test	0.021*	18:19:55
Air Blank	0.000	18:20:25
Air Blank	0.000	18:20:52
Subject Test	XXX**	18:21:57
Air Blank	0.000	18:22:31

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

Core Version: 8105.48

03/07/2013

18:22:58

5 of 7
TS4
SC

Standard Lot#: 23511100A2

Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847

TUCSON PD

Operator: T. GALLEGOS #28847

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	18:23:27
0.100 Cal Check	0.099	18:23:48
Air Blank	0.000	18:24:17

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

6077
TS4
SC
3/22/13

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 3-7-2013 TIME 1801
Intoxilyzer Serial # 80-000323 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

Core Version: 8105.48

03/07/2013

18:01:44

7 of 7
Bh

Standard Lot#: 23511100A2

Last Changed By: T. GALLEGOS #28847

SC

QAS: T. GALLEGOS #28847

TUCSON PD

Operator: T. GALLEGOS #28847

TUCSON PD

Test g/210L Time

Air Blank 0.000 18:02:14
0.100 Cal Check 0.100 18:02:35
Air Blank 0.000 18:03:04

Preventative Maintenance
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

SC
3/22/13

Date: 3-7-2013
 QAS: Terry Gallegos 28847
 Intoxilyzer Serial Number: 80-000323
 New Lot Number: 2351100A2
 New Expiration Date: 11-1-2013

10/2
TS4

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
 Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
 Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove ~~and sign~~ printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

Core Version: 8105.48

03/07/2013

17:59:01

2 of 2
TSh

50

Standard Lot#: 23511100A2

Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847

TUCSON PD

Operator: T. GALLEGOS #28847

TUCSON PD

Test g/210L Time

Air Blank 0.000 17:59:31
0.100 Cal Check 0.100 17:59:52
Air Blank 0.000 18:00:22

Preventative Maintenance

Performed -- Timer Reset #28847

QAS: T. GALLEGOS #28847