



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
2.10.2015

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/4/2015 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 1/8/2015 and 2/4/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on \_\_\_\_\_.  
The instrument was returned to service on \_\_\_\_\_.

William L Loomis

William L. Loomis


2/5/2015

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 5<sup>th</sup>, day of February, 2015.

Elma Escarcega

Notary Public

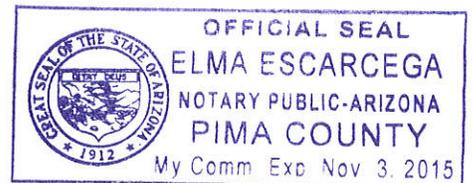

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
WLL  
TSY  
2-10-2015

QA SPECIALIST WILL LOOMIS 53545 AGENCY Tucson Police Department

DATE 2/4/2015 TIME 16:02

INTOXILYZER SERIAL # 80-000323 LOCATION TOXICOLOGY - CRIME LAB

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Will Loomis

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
Core Version: 8105.48  
02/04/2015 16:04:07

2 of 7  
WLL  
154

Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1980  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:05:00
Diagnostic Test	Pass	16:05:34
Air Blank	0.000	16:06:02
0.100 Cal Check	0.099	16:06:23
Air Blank	0.000	16:06:52
Subject Test	0.000	16:07:08
Air Blank	0.000	16:07:36
Wait	ABT*	
Air Blank	0.000	16:08:11

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
Core Version: 8105.48  
02/04/2015

16:08:20

3 of 7  
WLL  
B4

Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1980  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:09:11
Diagnostic Test	Pass	16:09:45
Air Blank	0.000	16:10:13
0.100 Cal Check	RFI*	16:10:28
Air Blank	IPS**	16:10:37

\*RFI Detect

\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
Core Version: 8105.48  
02/04/2015

16:10:51

4 of 7  
WLL  
TS4

---

Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: W. LOOMIS #53545  
TUCSON PD

---

Operator: W. LOOMIS #53545  
TUCSON PD

---

Subject: SQAP,TEST

DOB: 11/11/1980

Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	16:11:48
Diagnostic Test	Pass	16:12:22
Air Blank	0.000	16:12:49
0.100 Cal Check	0.099	16:13:10
Air Blank	0.000	16:13:40
Subject Test	0.486*	16:16:47
Air Blank	0.000	16:17:21
Air Blank	0.000	16:17:48
Subject Test	XXX**	16:18:12
Air Blank	0.000	16:18:46

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
Core Version: 8105.48  
02/04/2015 16:19:15

5 of 7  
WLL  
BY

-----  
Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	16:19:45
0.100 Cal Check	0.098	16:20:06
Air Blank	0.000	16:20:35

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
WLL  
TS4  
2.10.2015

QA Specialist WILL LOOMIS 53945 AGENCY Tucson Police Department  
DATE 2/4/2015 TIME 16:02  
Intoxilyzer Serial # 80-000323 Location TOXICOLOGY - CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution 11 AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
Core Version: 8105.48  
02/04/2015 16:02:08

7 of 7  
WLL  
154

-----  
Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	16:02:38
0.100 Cal Check	0.099	16:02:59
Air Blank	0.000	16:03:28

Preventative Maintenance  
Performed -- Timer Reset