



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TG  
7-6-2015

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/28/2015 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 4/30/2015 and 5/28/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.

William L Loomis

William L. Loomis


6/8/2015

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 8th, day of June, 20 15.

Elma Escarcega

Notary Public

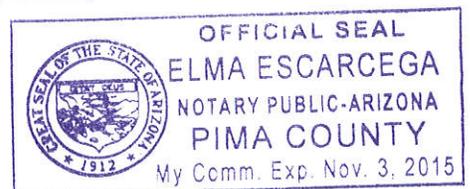

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
WLL  
TS4 7-6-2015

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 5/28/2015 TIME 14:26

INTOXILYZER SERIAL # 80-000323 LOCATION CRIME LAB

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Will Loomis

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
Core Version: 8105.48  
05/28/2015 14:28:58

Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP,TEST  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:29:48
Diagnostic Test	Pass	14:30:23
Air Blank	0.000	14:30:50
0.100 Cal Check	0.099	14:31:11
Air Blank	0.000	14:31:40
Subject Test	0.000	14:31:58
Air Blank	0.000	14:32:25
Wait	ABT*	
Air Blank	0.000	14:32:59

\*Sequence Aborted

2 of 7  
wll  
TS4

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
Core Version: 8105.48  
05/28/2015 14:33:10

3 of 7  
WU  
T84

---

Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: W. LOOMIS #53545  
TUCSON PD

---

Operator: W. LOOMIS #53545  
TUCSON PD

---

Subject: SQAP,TEST  
DOB: 11/11/1987  
Sex: M Weight: 222

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:34:05
Diagnostic Test	Pass	14:34:39
Air Blank	0.000	14:35:07
0.100 Cal Check	RFI*	14:35:22
Air Blank	IPS**	14:35:29

\*RFI Detect  
\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
Core Version: 8105.48  
05/28/2015 14:35:37

4 of 7  
WLV  
TS4

---

Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: W. LOOMIS #53545  
TUCSON PD

---

Operator: W. LOOMIS #53545  
TUCSON PD

---

Subject: SQAP,TEST  
DOB: 11/11/1987  
Sex: M Weight: 222

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:36:26
Diagnostic Test	Pass	14:37:00
Air Blank	0.000	14:37:27
0.100 Cal Check	0.098	14:37:48
Air Blank	0.000	14:38:18
Subject Test	0.268*	14:41:25
Air Blank	0.000	14:41:58
Air Blank	0.000	14:42:26
Subject Test	XXX**	14:42:53
Air Blank	0.000	14:43:27

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
Core Version: 8105.48  
05/28/2015 14:43:48

-----  
Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:44:17
0.100 Cal Check	0.097	14:44:37
Air Blank	0.000	14:45:07

Preventative Maintenance  
Performed -- Timer Reset

5 of 7  
TJH  
WU

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
WLL  
TSY 7.6.2015

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department  
DATE 5/28/2015 TIME 14:26  
Intoxilyzer Serial # 80-000323 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
Core Version: 8105.48  
05/28/2015 14:26:33

-----  
Standard Lot#: 25913100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:27:02
0.100 Cal Check	0.098	14:27:23
Air Blank	0.000	14:27:53

Preventative Maintenance  
Performed -- Timer Reset

7 of 7  
WR  
TS4