



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
2-11-2016

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/4/2016 Cal 1 0.102 AC Cal 2 0.101 AC

Between the dates of 1/7/2016 and 2/4/2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on \_\_\_\_\_.  
The instrument was returned to service on \_\_\_\_\_.

William L Loomis

William L. Loomis

2/10/2016

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 10<sup>th</sup> day of February, 2016.

Elma Escarcega

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

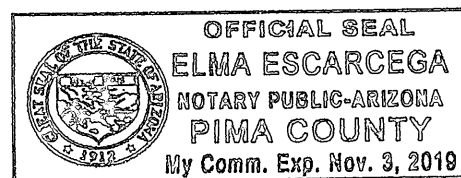


EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
wlv  
By  
2.11.2016

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 2/4/2016 TIME 13:51

INTOXILYZER SERIAL # 80-000323 LOCATION CRIME LAB

( ☒ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ☒ ) 1. Clock time check.  
( ☒ ) 2. Date check.

OPERATIONAL TESTS

- ( ☒ ) 1. Alcohol-free subject test result 0. 000 AC.  
( ☒ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
( ☒ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
( ☒ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 101 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Willie Loomis

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

Core Version: 8105.48

02/04/2016

13:54:00

Standard Lot#: 12415100A4

Last Changed By: W. LOOMIS

#53545

QAS: W. LOOMIS

#53545

TUCSON PD

Operator: W. LOOMIS

#53545

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1987

Sex: M

Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:54:51
Diagnostic Test	Pass	13:55:25
Air Blank	0.000	13:55:52
0.100 Cal Check	0.102	13:56:13
Air Blank	0.000	13:56:42
Subject Test	0.000	13:57:00
Air Blank	0.000	13:57:27
Wait	ABT*	
Air Blank	0.000	13:58:02

\*Sequence Aborted

2 of 7  
WLL  
TS4

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

Core Version: 8105.48

02/04/2016

13:58:08

3 of 7  
wu  
B4

Standard Lot#: 12415100A4

Last Changed By: W. LOOMIS

#53545

QAS: W. LOOMIS

#53545

TUCSON PD

Operator: W. LOOMIS

#53545

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1987

Sex: M

Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:58:58
Diagnostic Test	Pass	13:59:32
Air Blank	0.000	13:59:59
0.100 Cal Check	RFI*	14:00:14
Air Blank	IPS**	14:00:22

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

Core Version: 8105.48

02/04/2016

14:00:56

4 of 7  
ms  
TS4

Standard Lot#: 12415100A4

Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545

TUCSON PD

Operator: W. LOOMIS #53545

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1987

Sex: M

Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:01:48
Diagnostic Test	Pass	14:02:22
Air Blank	0.000	14:02:49
0.100 Cal Check	0.102	14:03:10
Air Blank	0.000	14:03:39
Subject Test	0.183*	14:06:47
Air Blank	0.000	14:07:20
Air Blank	0.000	14:07:47
Subject Test	XXX**	14:08:04
Air Blank	0.000	14:08:38

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

Core Version: 8105.48

02/04/2016

14:08:51

-----  
Standard Lot#: 12415100A4

Last Changed By: W. LOOMIS #53545  
-----

QAS: W. LOOMIS #53545

TUCSON PD  
-----

Operator: W. LOOMIS #53545

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:09:20
0.100 Cal Check	0.101	14:09:41
Air Blank	0.000	14:10:10

Preventative Maintenance  
Performed -- Timer Reset

5 of 7  
WV  
By

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
WLV  
104  
2-11-2016

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department  
DATE 2/4/2016 TIME 13:57  
Intoxilyzer Serial # 80-000323 Location CRIME LAB

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 102 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Willie Loomis

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

Core Version: 8105.48

02/04/2016

13:51:46

-----  
Standard Lot#: 12415100A4

Last Changed By: W. LOOMIS #53545  
-----

QAS: W. LOOMIS #53545

TUCSON PD  
-----

Operator: W. LOOMIS #53545

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:52:15
0.100 Cal Check	0.102	13:52:36
Air Blank	0.000	13:53:05

Preventative Maintenance

Performed -- Timer Reset

7 of 7  
WU  
TS4