



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

B4
8.22.2016

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8/18/2016 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 7/21/2016 and 8/18/2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

William L Loomis

William L. Loomis


8/22/2016

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 22nd, day of August, 2016.

Elma Escarcega

Notary Public

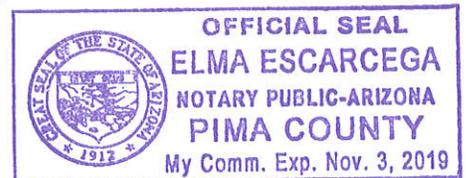

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

*1 of 8 will
by 8.22.2016*

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 8/18/2016 TIME 12:19

INTOXILYZER SERIAL # 80-000323 LOCATION CRIME LAB

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Willie Loomis

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
08/18/2016

12:22:36

Standard Lot#: 12415100A4
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 12:23:30 |
| Diagnostic Test | Fail* | 12:23:35 |
| Air Blank | 0.000 | 12:24:03 |

* Voltage/Current Test

2 of 8
WLL
TS4

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
08/18/2016

12:25:24

Standard Lot#: 12415100A4
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 12:26:44 |
| Diagnostic Test | Pass | 12:27:19 |
| Air Blank | 0.000 | 12:27:46 |
| 0.100 Cal Check | 0.099 | 12:28:07 |
| Air Blank | 0.000 | 12:28:36 |
| Subject Test | 0.000 | 12:28:53 |
| Air Blank | 0.000 | 12:29:21 |
| Wait | ABT* | |
| Air Blank | 0.000 | 12:29:57 |

*Sequence Aborted

3 of 8
WU
T84

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
08/18/2016

12:30:11

Standard Lot#: 12415100A4
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987

Sex: M Weight: 222

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 12:31:32 |
| Diagnostic Test | Pass | 12:32:06 |
| Air Blank | 0.000 | 12:32:34 |
| 0.100 Cal Check | 0.099 | 12:32:55 |
| Air Blank | RFI* | 12:33:01 |
| Air Blank | IPS** | 12:33:10 |

*RFI Detect

**Improper Sample

4 of 8
WLL
By

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
08/18/2016

12:33:36

Standard Lot#: 12415100A4
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 12:34:56 |
| Diagnostic Test | Pass | 12:35:30 |
| Air Blank | 0.000 | 12:35:57 |
| 0.100 Cal Check | 0.100 | 12:36:18 |
| Air Blank | 0.000 | 12:36:47 |
| Subject Test | 0.094* | 12:39:55 |
| Air Blank | 0.000 | 12:40:27 |
| Air Blank | 0.000 | 12:40:54 |
| Subject Test | XXX** | 12:41:13 |
| Air Blank | 0.000 | 12:41:46 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

5 of 8
WLL
Bh

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
Core Version: 8105.48
08/18/2016 12:42:01

Standard Lot#: 12415100A4
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 12:42:30 |
| 0.100 Cal Check | 0.099 | 12:42:51 |
| Air Blank | 0.000 | 12:43:21 |

Preventative Maintenance
Performed -- Timer Reset

6 of 8
TS4
WLL

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

7 of 8
WLL
By 8.22.2016

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department
DATE 8/18/2016 TIME 12:19
Intoxilyzer Serial # 80-000323 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE William Loomis

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

