



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
9.22.2016

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/15/2016 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 8/18/2016 and 9/15/2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on \_\_\_\_\_.  
The instrument was returned to service on \_\_\_\_\_.

William L Loomis

William L. Loomis


9/19/2016

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 19th, day of September, 2016.

Elma Escarcega

Notary Public

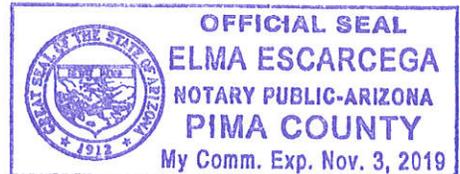

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
WLL  
By 9.22.2016

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 9/15/2016 TIME 13:14

INTOXILYZER SERIAL # 80-000323 LOCATION CRIME LAB

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Willie Loomis

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
Core Version: 8105.48  
09/15/2016

13:16:32

2 of 7  
WLL  
TBh

Standard Lot#: 12415100A4  
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:17:25
Diagnostic Test	Pass	13:17:59
Air Blank	0.000	13:18:26
0.100 Cal Check	0.099	13:18:47
Air Blank	0.000	13:19:16
Subject Test	0.000	13:19:35
Air Blank	0.000	13:20:02
Wait	ABT*	
Air Blank	0.000	13:20:37

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
Core Version: 8105.48  
09/15/2016 13:20:44

Standard Lot#: 12415100A4  
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:21:46
Diagnostic Test	Pass	13:22:20
Air Blank	0.000	13:22:47
0.100 Cal Check	RFI*	13:23:01
Air Blank	IPS**	13:23:09

\*RFI Detect  
\*\*Improper Sample

3 of 7  
WLL  
B4

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
Core Version: 8105.48  
09/15/2016 13:23:32

Standard Lot#: 12415100A4  
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545  
TUCSON PD

Operator: W. LOOMIS #53545  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:24:22
Diagnostic Test	Pass	13:24:56
Air Blank	0.000	13:25:24
0.100 Cal Check	0.099	13:25:45
Air Blank	0.000	13:26:14
Subject Test	0.033*	13:29:21
Air Blank	0.000	13:29:53
Air Blank	0.000	13:30:21
Subject Test	XXX**	13:30:34
Air Blank	0.000	13:31:08

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

4 of 7  
WU  
TS4

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
Core Version: 8105.48  
09/15/2016 13:31:22

-----  
Standard Lot#: 12415100A4  
Last Changed By: W. LOOMIS #53545  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:31:52
0.100 Cal Check	0.099	13:32:13
Air Blank	0.000	13:32:43

Preventative Maintenance  
Performed -- Timer Reset

5 of 7  
WLL  
BSG

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
WLL  
TSL 9.22.2016

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department  
DATE 9/15/2016 TIME 13:14  
Intoxilyzer Serial # 80-000323 Location CRIME LAB

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
Core Version: 8105.48  
09/15/2016 13:14:17

-----  
Standard Lot#: 12415100A4  
Last Changed By: W. LOOMIS #53545  
-----

QAS: W. LOOMIS #53545  
TUCSON PD  
-----

Operator: W. LOOMIS #53545  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:14:46
0.100 Cal Check	0.099	13:15:07
Air Blank	0.000	13:15:36

Preventative Maintenance  
Performed -- Timer Reset

7 of 7  
DWL  
TBH