



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

ju

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-22-2006 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 12-1-2006 and 12-22-2006 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

Terry S. Gallegos


Terry S. Gallegos

1-4-2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 4, day of

January, 2007.

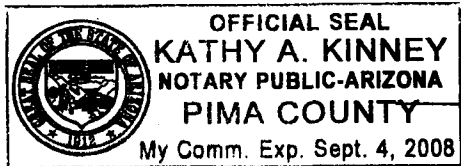

Kathy A. Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28844 AGENCY Tucson Police DepartmentDATE 12.22.2006 TIME 1534INTOXILYZER SERIAL # 80-000331 LOCATION ODW

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
-
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.
- 000
- AC.
-
- (✓) 2. Error recognition logic system functioning. ✓
-
- Not a Successfully Completed Test Sequence printed
-
- (✓) 3. Proper sample recognition system.
-
- Not a Successfully Completed Test Sequence printed
-
- Deficient sample printed.
-
- (✓) 4. Standard Calibration Check standard 0.
- 100
- AC
-
- Result: 0.
- 100
- AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Terry Gallegos

2058
TS4

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000331

12/22/2006

15:37:12

Standard Lot#: 502602I

Last Changed By: J. COOMBS

#46160

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:38:26
Diagnostic Test	Pass	15:38:54
Air Blank	0.000	15:39:21
0.100 Cal Check	0.100	15:39:40
Air Blank	0.000	15:40:10
Subject Test	0.000	15:40:33
Air Blank	0.000	15:41:00
Five Minute Wait	ABT*	
Air Blank	0.000	15:41:51

*Sequence Aborted

Not a Successfully Completed Test Sequence

3058
TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
12/22/2006

15:44:44

Standard Lot#: 502602I
Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	15:45:25
Air Blank	0.000	15:45:54

*Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000331

12/22/2006

15:46:12

Standard Lot#: 502602I

Last Changed By: J. COOMBS

#46160

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	15:46:45
Air Blank	0.000	15:47:13

*RFI Detect

Not a Successfully Completed Test Sequence

5080
TSh

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000331

12/22/2006

15:47:52

Standard Lot#: 502602I

Last Changed By: J. COOMBS

#46160

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:48:47
Diagnostic Test	Pass	15:49:14
Air Blank	0.000	15:49:41
0.100 Cal Check	0.100	15:50:00
Air Blank	0.000	15:50:30
Subject Test	0.023*	15:53:38
Air Blank	0.000	15:54:09
Air Blank	0.000	15:54:36
Subject Test	XXX**	15:55:12
Air Blank	0.000	15:55:45

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000331

12/22/2006

15:57:02

Standard Lot#: 502602I

Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	15:57:30
Cal Check	0.100	15:57:48
Air Blank	0.000	15:58:18


Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 12-22-2006 TIME 1534
Intoxilyzer Serial # 80-000331 Location ODW

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE *Terry Gallegos*

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.