



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

64

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/19/07 Cal 1 0.102 AC Cal 2 0.101 AC

Between the dates of 12/22/06 and 1/19/07 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Jennifer L. Coombs  
Jennifer L. Coombs

2/9/07

ADPS QAS Permit # 33120

SUBSCRIBED AND SWORN to before me this 9, day of February, 2007.

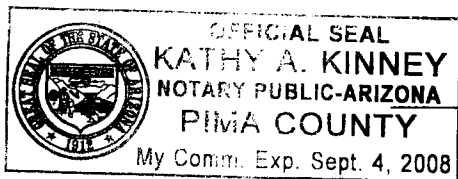

Kathy A. Kinney  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Jennifer Coombs 46160 AGENCY Tucson Police Department  
DATE January 19, 2007 TIME 1751  
INTOXILYZER SERIAL # 80-000331 LOCATION ODW

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.  
(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0.000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0.100 AC  
Result: 0.101 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS Gas tank change prior to G-2.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Jennifer L. Coombs

DPS Form Exh G-4 (Rev. 05-01)

INTOXILYZER 8000 .

Location: TUCSON POLICE DEPT

Serial Number: 80-000331

01/19/2007

18:11:00

-----  
Standard Lot#: 629301I

Last Changed By: J. COOMBS #46160  
-----

QAS: J. COOMBS #46160

TUCSON PD CRIME LAB  
-----

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | 0.000  | 18:11:27 |
| Cal Check | 0.101  | 18:11:46 |
| Air Blank | 0.000  | 18:12:17 |

  
\_\_\_\_\_  
Operator's Signature

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000331

01/19/2007

17:53:58

Standard Lot#: 629301I

Last Changed By: J. COOMBS

#46160

QAS: J. COOMBS

#46160

TUCSON PD CRIME LAB

Operator: J. COOMBS

#46160

TUCSON PD CRIME LAB

Subject: SQAP,TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

| Test             | g/210L | Time     |
|------------------|--------|----------|
| Air Blank        | 0.000  | 17:54:54 |
| Diagnostic Test  | Pass   | 17:55:21 |
| Air Blank        | 0.000  | 17:55:48 |
| 0.100 Cal Check  | 0.102  | 17:56:07 |
| Air Blank        | 0.000  | 17:56:38 |
| Subject Test     | 0.000  | 17:57:18 |
| Air Blank        | 0.000  | 17:57:46 |
| Five Minute Wait | ABT*   |          |
| Air Blank        | 0.000  | 17:58:25 |

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000331  
01/19/2007

17:58:44

---

Standard Lot#: 629301I  
Last Changed By: J. COOMBS #46160

---

QAS: J. COOMBS #46160  
TUCSON PD CRIME LAB

---

Operator: J. COOMBS #46160  
TUCSON PD CRIME LAB

---

Subject: SQAP, TEST

DOB: 11/11/1111

Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | IPS*   | 17:59:22 |
| Air Blank | 0.000  | 17:59:50 |

\*Improper Sample

**Not a Successfully Completed Test Sequence**

*Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.*

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000331

01/19/2007

17:59:57

Standard Lot#: 629301I

Last Changed By: J. COOMBS

#46160

QAS: J. COOMBS

#46160

TUCSON PD CRIME LAB

Operator: J. COOMBS

#46160

TUCSON PD CRIME LAB

Subject: SOAP,TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 18:00:52 |
| Diagnostic Test | Pass   | 18:01:19 |
| Air Blank       | 0.000  | 18:01:46 |
| 0.100 Cal Check | RFI*   | 18:02:02 |
| Air Blank       | 0.000  | 18:02:33 |

\*RFI Detect

**Not a Successfully Completed Test Sequence**

*Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department*

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000331

01/19/2007

18:02:40

Standard Lot#: 629301I

Last Changed By: J. COOMBS

#46160

QAS: J. COOMBS

#46160

TUCSON PD CRIME LAB

Operator: J. COOMBS

#46160

TUCSON PD CRIME LAB

Subject: TEST,TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 18:03:38 |
| Diagnostic Test | Pass   | 18:04:06 |
| Air Blank       | 0.000  | 18:04:33 |
| 0.100 Cal Check | 0.101  | 18:04:52 |
| Air Blank       | 0.000  | 18:05:23 |
| Subject Test    | 0.027* | 18:08:31 |
| Air Blank       | 0.000  | 18:09:02 |
| Air Blank       | 0.000  | 18:09:29 |
| Subject Test    | XXX**  | 18:10:02 |
| Air Blank       | 0.000  | 18:10:34 |

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

## STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Jennifer Coombs 46160 AGENCY Tucson Police Department  
DATE January 19, 2007 TIME 1751  
Intoxilyzer Serial # 80-000331 Location ODW

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Jennifer L. Coombs

DPS Form Exh G-2 (Rev. 05-01)



INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000331

01/19/2007

17:51:05


-----  
Standard Lot#: 629301I

Last Changed By: J. COOMBS #46160  
-----

QAS: J. COOMBS #46160

TUCSON PD CRIME LAB  
-----

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | 0.000  | 17:51:33 |
| Cal Check | 0.102  | 17:51:51 |
| Air Blank | 0.000  | 17:52:22 |

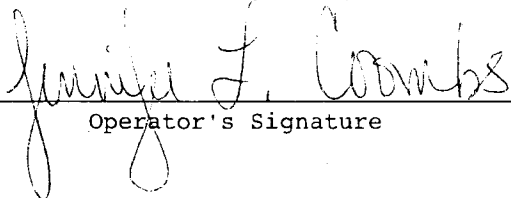  
\_\_\_\_\_  
Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

1 of 1  
jlc

INTOXILYZER SERVICE RECORD <sup>TS4</sup>  
FOR  
CALIBRATION GAS TANK CHANGE

Date: 1/19/07  
QAS: Jennifer Coombs 46160  
Intoxilyzer Serial Number: 80-000331  
New Lot Number: 629301 I  
New Expiration Date: 10/23/2008

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
Enter Target Value: 0.100  
Enter Control Lot #: From the tank label  
Enter Exp Date: From the tank label (mm/dd/yyyy)  
Edit info Y/N? Y or N  
Update Std Info? Y ENTER  
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 10.2 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE

Jennifer L. Coombs