



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

156

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/15/07 Cal 1 0.101 AC Cal 2 0.101 AC

Between the dates of 1/19/07 and 2/15/07 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Jennifer L. Coombs
Jennifer L. Coombs

2/16/07

ADPS QAS Permit # 33120

SUBSCRIBED AND SWORN to before me this 16, day of February, 2007.

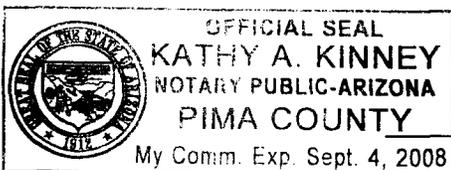

Kathy A. Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

1 of 8
file

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

BH

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Jennifer Coombs 46160 AGENCY Tucson Police Department

DATE February 15, 2007 TIME 1040

INTOXILYZER SERIAL # 80-000331 LOCATION ODW

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed ✓✓
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Jennifer L. Coombs

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
02/15/2007

10:54:54

BH

2 of 8
file

Standard Lot#: 629301I
Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160
TUCSON PD CRIME LAB

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:55:47
Diagnostic Test	Pass	10:56:15
Air Blank	IPS*	10:56:23
Air Blank	0.000	10:56:51

*Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3078
fle

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
02/15/2007 10:42:15

B4

Standard Lot#: 629301I
Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160
TUCSON PD CRIME LAB

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	ABT*	10:43:15
Air Blank	0.000	10:43:43

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
02/15/2007

10:44:52

4 of 8
JL

Standard Lot#: 629301I
Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160
TUCSON PD CRIME LAB

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:45:51
Diagnostic Test	Pass	10:46:19
Air Blank	0.000	10:46:46
0.100 Cal Check	0.102	10:47:05
Air Blank	0.000	10:47:34
Subject Test	0.000	10:47:50
Air Blank	0.000	10:48:17
Five Minute Wait		
Air Blank	0.000	10:53:17
Subject Test	RFI*	10:53:55
Air Blank	0.000	10:54:23

*RFI Detect

Not a Successfully Completed Test Sequence

5 of 8
file

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
02/15/2007 10:56:55

B4

Standard Lot#: 629301I
Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160
TUCSON PD CRIME LAB

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:57:58
Diagnostic Test	Pass	10:58:25
Air Blank	0.000	10:58:52
0.100 Cal Check	0.101	10:59:11
Air Blank	0.000	10:59:42
Subject Test	0.156*	11:02:50
Air Blank	0.000	11:03:22
Air Blank	0.000	11:03:49
Subject Test	XXX**	11:04:05
Air Blank	0.000	11:04:38

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
02/15/2007 11:05:50


Loofg
JL

Standard Lot#: 629301I
Last Changed By: J. COOMBS #46160

BA

QAS: J. COOMBS #46160
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	11:06:18
Cal Check	0.101	11:06:37
Air Blank	0.000	11:07:07

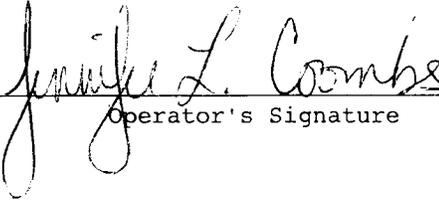

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

*7 of 8
file*

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

164

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Jennifer Coombs 46160 AGENCY Tucson Police Department
DATE February 15, 2007 TIME 1040
Intoxilyzer Serial # 8000031 Location ODW

- () 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- () 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- () 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- () 4. Air blank completed.
- () 5. Calibration check completed. Test results 0. 101 AC.
- () 6. Air blank completed.
- () 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Jennifer L. Coombs

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
02/15/2007 10:40:52

*Soft
file*

B4

Standard Lot#: 629301I
Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	10:41:19
Cal Check	0.101	10:41:38
Air Blank	0.000	10:42:08

Jennifer L. Coombs

Operator's Signature