



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

Bh
7-10-2007

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7-5-07 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 6-7-07 and 7-5-07 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

7-9-07

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 9, day of July, 2007.

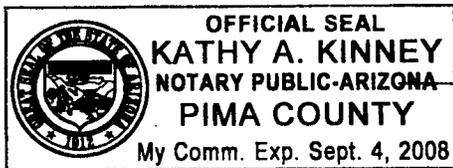

Kathy A Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

154

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 7-5-07 TIME 0828

INTOXILYZER SERIAL # 80-000331 LOCATION DDW

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
07/05/2007

08:29:39

Ph

Standard Lot#: 629301I
Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:30:34
Diagnostic Test	Pass	08:31:01
Air Blank	0.000	08:31:28
0.100 Cal Check	0.100	08:31:47
Air Blank	0.000	08:32:18
Subject Test	0.000	08:32:35
Air Blank	0.000	08:33:02
Five Minute Wait	ABT*	
Air Blank	0.000	08:33:36

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 9 cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
07/05/2007

08:34:23

Ph

Standard Lot#: 629301I
Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:35:17
Diagnostic Test	Pass	08:35:44
Air Blank	0.000	08:36:11
0.100 Cal Check	RFI*	08:36:24
Air Blank	0.000	08:36:55

*RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 9 cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
07/05/2007 08:37:03

T84

Standard Lot#: 629301I
Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	08:37:40
Air Blank	0.000	08:38:08

*Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5069 cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
07/05/2007

08:45:46

B4

Standard Lot#: 629301I
Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:46:42
Diagnostic Test	Pass	08:47:09
Air Blank	0.000	08:47:37
0.100 Cal Check	0.100	08:47:55
Air Blank	0.000	08:48:26
Subject Test	0.023*	08:51:34
Air Blank	0.000	08:52:05
Air Blank	0.000	08:52:33
Subject Test	ABT**	08:52:43
Air Blank	0.000	08:53:11

*Deficient Sample
**Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

6 of 9 cu

Tsh

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
07/05/2007 08:40:01

Standard Lot#: 629301I
Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:40:53
Diagnostic Test	Pass	08:41:20
Air Blank	0.000	08:41:48
0.100 Cal Check	0.100	08:42:06
Air Blank	0.000	08:42:38
Subject Test	XXX*	08:44:24
Air Blank	0.000	08:44:58

*Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

749 cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
07/05/2007 08:53:49

TS4

Standard Lot#: 629301I
Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:54:16
Cal Check	0.100	08:54:35
Air Blank	0.000	08:55:06

Claire Conley

Operator's Signature

8 of 9 cu

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 7-5-07 TIME 0828
Intoxilyzer Serial # 80-000331 Location ODW

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

989 cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
07/05/2007 08:28:15

TS4

Standard Lot#: 629301I
Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:28:42
Cal Check	0.100	08:29:01
Air Blank	0.000	08:29:32

Clare Conley

Operator's Signature