



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

TS4  
8-8-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8-1-07 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 7-5-07 and 8-1-07 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Clare Conley


Clare L. Conley

8-3-07

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 3, day of

August, 2007.

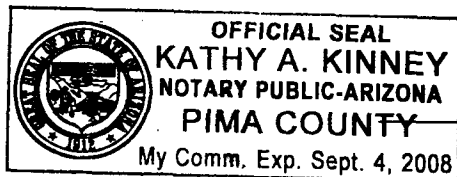

Kathy A Kinney

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

B4  
8-8-07

QA SPECIALIST Clare Conley AGENCY Tucson Police Department  
DATE 8-1-07 TIME 10:37  
INTOXILYZER SERIAL # 80-000331 LOCATION ODW

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.  
( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.  
( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ✓ NO     

COMMENTS None.

SIGNATURE Clare Conley

2 of 8 m

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000331  
08/01/2007

10:38:58

154

Standard Lot#: 629301I  
Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:39:49
Diagnostic Test	Pass	10:40:16
Air Blank	0.000	10:40:44
0.100 Cal Check	0.100	10:41:03
Air Blank	0.000	10:41:33
Subject Test	0.000	10:41:48
Air Blank	ABT*	10:41:55
Air Blank	0.000	10:42:23

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 8 Cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000331

08/01/2007

10:42:34

BS4

Standard Lot#: 629301I

Last Changed By: J. COOMBS

#46160

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:43:19
Air Blank	0.000	10:43:48

\*RFI Detect

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 8 Cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000331

08/01/2007

10:43:52

154

Standard Lot#: 629301I

Last Changed By: J. COOMBS

#46160

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	10:44:26
Air Blank	0.000	10:44:54

\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

548Cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000331

08/01/2007

10:44:58

104

Standard Lot#: 629301I

Last Changed By: J. COOMBS

#46160

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:45:51
Diagnostic Test	Pass	10:46:18
Air Blank	0.000	10:46:46
0.100 Cal Check	0.100	10:47:04
Air Blank	0.000	10:47:35
Subject Test	0.362*	10:50:44
Air Blank	0.000	10:51:16
Air Blank	0.000	10:51:43
Subject Test	XXX**	10:52:02
Air Blank	0.000	10:52:36

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

6 of 8 cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000331

08/01/2007

10:52:55

B4

-----  
Standard Lot#: 629301I

Last Changed By: J. COOMBS #46160  
-----

QAS: C. CONLEY #34517

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:53:23
Cal Check	0.099	10:53:41
Air Blank	0.000	10:54:12

*C. Conley*

\_\_\_\_\_  
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 8-1-07 TIME 10:37  
Intoxilyzer Serial # 80-000331 Location ODW

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 100 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Clare Conley

DPS Form Exh G-2 (Rev. 05-01)



8 of 8 cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000331

08/01/2007

10:37:36

BH

-----  
Standard Lot#: 629301I

Last Changed By: J. COOMBS #46160  
-----

QAS: C. CONLEY #34517

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:38:03
Cal Check	0.100	10:38:22
Air Blank	0.000	10:38:53

*Clare Conley*

\_\_\_\_\_  
Operator's Signature