



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

TSy  
9.4.2007

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8-30-07 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 8-1-07 and 8-30-07 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

8.31.07

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 31, day of August, 2007.

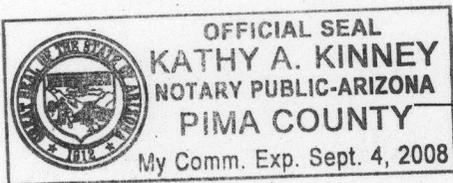

Kathy A Kinney

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

1069 cm  
154

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 8-30-07 TIME 10:07

INTOXILYZER SERIAL # 80-000331 LOCATION ODW

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 100 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000331  
08/30/2007

10:11:36

2069 cu  
Bh

Standard Lot#: 719801I  
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:12:27
Diagnostic Test	Pass	10:12:55
Air Blank	0.000	10:13:22
0.100 Cal Check	0.101	10:13:41
Air Blank	0.000	10:14:11
Subject Test	0.000	10:14:28
Air Blank	ABT*	10:14:34
Air Blank	0.000	10:15:02

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000331  
08/30/2007 10:15:06

389a  
by

---

Standard Lot#: 719801I  
Last Changed By: C. CONLEY #34517

---

QAS: C. CONLEY #34517  
TUCSON PD

---

Operator: C. CONLEY #34517  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	10:15:41
Air Blank	IPS**	10:15:55

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000331  
08/30/2007

10:16:01

459 cu  
B4

---

Standard Lot#: 719801I  
Last Changed By: C. CONLEY #34517

---

QAS: C. CONLEY #34517  
TUCSON PD

---

Operator: C. CONLEY #34517  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:16:48
Diagnostic Test	Pass	10:17:16
Air Blank	0.000	10:17:43
0.100 Cal Check	0.101	10:18:02
Air Blank	0.000	10:18:32
Subject Test	0.165*	10:21:40
Air Blank	0.000	10:22:12
Air Blank	0.000	10:22:40
Subject Test	XXX**	10:23:05
Air Blank	0.000	10:23:39

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000331  
08/30/2007 10:23:58

5 of 9 cu  
BH

-----  
Standard Lot#: 719801I  
Last Changed By: C. CONLEY #34517  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:24:25
Cal Check	0.100	10:24:44
Air Blank	0.000	10:25:14

*Clare Conley*

\_\_\_\_\_  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6-269 cu  
104

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 8-30-07 TIME 10:07  
Intoxilyzer Serial # 80-000331 Location DDW

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0.100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0.101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000331  
08/30/2007 10:09:12

7069 cu  
TS4

-----  
Standard Lot#: 719801I  
Last Changed By: C. CONLEY #34517  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:09:40
Cal Check	0.101	10:09:58
Air Blank	0.000	10:10:29

*Clare Conley*

Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

8569 cu  
154

Date: 8-30-07  
 QAS: Clare Conley  
 Intoxilyzer Serial Number: 80-000331  
 New Lot Number: 7198011  
 New Expiration Date: 23 Jul 09

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change." Select "D" for dry gas and enter the following information:  
 Enter Target Value: 0.100  
 Enter Control Lot #: From the tank label  
 Enter Exp Date: From the tank label (mm/dd/yyyy)  
 Edit info Y/N? Y or N  
 Update Std Info? Y ENTER  
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing." Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 101 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000331  
08/30/2007 10:07:35

9 of 9 cu  
PH

-----  
Standard Lot#: 719801I  
Last Changed By: C. CONLEY #34517  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:08:02
Cal Check	0.101	10:08:21
Air Blank	0.000	10:08:51

*Clare Conley*

-----  
Operator's Signature