



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

BH
6.6.08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-5-08 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 5-8-08 and 6-5-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A
The instrument was returned to service on N/A


Clare Conley

Clare L. Conley

6.6.08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 6, day of June, 2008.

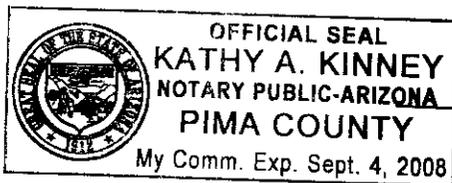

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7
Cn
By
6-6-08

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department
DATE 6-5-08 TIME 0905
INTOXILYZER SERIAL # 80-050331 LOCATION DUI HQ Team 2

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7 cu

BM

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
06/05/2008

09:08:32

Standard Lot#: 719801I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:09:26
Diagnostic Test	Pass	09:09:54
Air Blank	0.000	09:10:21
0.100 Cal Check	0.100	09:10:40
Air Blank	0.000	09:11:09
Subject Test	0.000	09:11:26
Air Blank	ABT*	09:11:36
Air Blank	0.000	09:12:04

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7 cu

ASU

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
06/05/2008 09:12:42

Standard Lot#: 719801I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:13:38
Diagnostic Test	Pass	09:14:05
Air Blank	RFI*	09:14:12
Air Blank	IPS**	09:14:26

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4.67 Cu

Tsh

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
06/05/2008

09:14:31

Standard Lot#: 719801I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:15:28
Diagnostic Test	Pass	09:15:56
Air Blank	0.000	09:16:23
0.100 Cal Check	0.099	09:16:42
Air Blank	0.000	09:17:12
Subject Test	0.080*	09:20:20
Air Blank	0.000	09:20:52
Air Blank	0.000	09:21:19
Subject Test	XXX**	09:21:42
Air Blank	0.000	09:22:15

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

5367 Cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
06/05/2008

09:22:41

TS4

Standard Lot#: 719801I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:23:08
Cal Check	0.099	09:23:27
Air Blank	0.000	09:23:57

Clare Conley

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6/5/7
cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

BA
6-6-08

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 6-5-08 TIME 0905
Intoxilyzer Serial # 80-00031 Location DUI HQ Team 2

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

7 of 7 Cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
06/05/2008 09:05:09

134

Standard Lot#: 719801I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:05:37
Cal Check	0.099	09:05:56
Air Blank	0.000	09:06:25

Clare Conley

QAS: C. Operator's Signature