



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TSM
12-24-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-19-08 Cal 1 0.102 AC Cal 20.100 AC

Between the dates of 11-20-08 and 12-19-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


12-23-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 23, day of December, 2008.

A.M. Hernandez

Notary Public

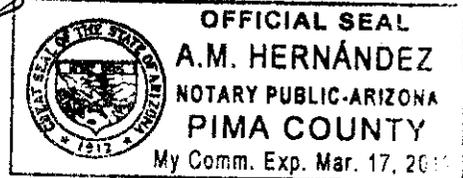

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1067
cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
12-24-08

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 12-19-08 TIME 0933

INTOXILYZER SERIAL # 80-000331 LOCATION Crime Laboratory

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
12/19/2008

09:35:15

2 of 7
Cu
TS

Standard Lot#: 811201I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:36:00
Diagnostic Test	Pass	09:36:28
Air Blank	0.000	09:36:55
0.100 Cal Check	0.101	09:37:14
Air Blank	0.000	09:37:44
Subject Test	0.000	09:38:07
Air Blank	ABT*	09:38:14
Air Blank	0.000	09:38:42

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
12/19/2008

09:39:01

3061

Cu

TSH

Standard Lot#: 811201I
Last Changed By: C. CONLEY- #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:39:32
Air Blank	IPS**	09:39:40

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
12/19/2008

09:39:50

487
Cu

Standard Lot#: 811201I
Last Changed By: C. CONLEY- #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

TS4

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:40:41
Diagnostic Test	Pass	09:41:08
Air Blank	0.000	09:41:36
0.100 Cal Check	0.102	09:41:54
Air Blank	0.000	09:42:25
Subject Test	0.249*	09:45:33
Air Blank	0.000	09:46:00
Air Blank	0.000	09:46:34
Subject Test	XXX**	09:46:54
Air Blank	0.000	09:47:28

*Deficient Sample

**Invalid Sample

Begin new deprivation period


Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
12/19/2008 09:47:55

Standard Lot#: 811201I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:48:22
Cal Check	0.100	09:48:41
Air Blank	0.000	09:49:11

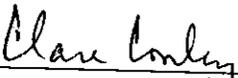

Operator's Signature

5067
Cu
TS4

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

b-67
Cu

754
12-24-08

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 12-19-08 TIME 0933
Intoxilyzer Serial # 80-000331 Location Crime Laboratory

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
12/19/2008 09:33:53

7567
CW

Standard Lot#: 811201I
Last Changed By: C. CONLEY #34517

154

QAS: C. CONLEY #34517 -
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:34:20
Cal Check	0.102	09:34:39
Air Blank	0.000	09:35:09

Clare Conley

Operator's Signature