



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 270 S. Stone Avenue Tucson, AZ 85701-1917

TSB
 3/17/09

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/13/09 Cal 1 0.104 AC Cal 2 0.102 AC

Between the dates of 2/12/09 and 3/13/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
 The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


3/17/09

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 17, day of March, 2009.

A.M. Hernandez

Notary Public

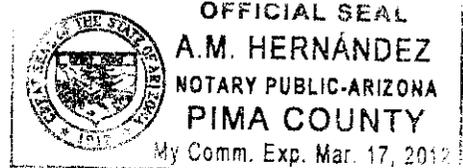

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

124
3-17-09

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 3/13/09 TIME 11:42

INTOXILYZER SERIAL # 80-000331 LOCATION Vehicle #432-Tm2

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 102 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7_{sc}

754

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
03/13/2009 11:43:39

Standard Lot#: 8112011
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:44:27
Diagnostic Test	Pass	11:44:56
Air Blank	0.000	11:45:23
0.100 Cal Check	0.100	11:45:42
Air Blank	0.000	11:46:12
Subject Test	0.000	11:46:27
Air Blank	0.000	11:46:54
Five Minute Wait	ABT*	
Air Blank	0.000	11:47:30

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

3 of 7 SC
154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
03/13/2009 11:47:44

Standard Lot#: 8112011
Last Changed By: C. CONLEY #34517

OAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	11:48:17
Air Blank	IPS**	11:48:26

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

4 of 7 SC

154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
03/13/2009 11:46:37

Standard Lot#: 8112011
Last Changed By: C. CONLEY #34517

GAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:49:30
Diagnostic Test	Pass	11:49:58
Air Blank	0.000	11:50:26
0.100 Cal Check	0.100	11:50:45
Air Blank	0.000	11:51:14
Subject Test	0.070*	11:54:23
Air Blank	0.000	11:54:54
Air Blank	0.000	11:55:21
Subject Test	XXX**	11:55:36
Air Blank	0.000	11:56:09

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully
Completed
Test Sequence

5 of 7 SC

154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
03/13/2009 11:56:26

Standard Lot#: 8112011
Last Changed By: C. CONLEY #34517

OAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:56:54
Cal Check	0.102	11:57:12
Air Blank	0.000	11:57:43

Sarah Cione
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 SC

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

TS4
3-17-09

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 3/13/09 TIME 11:42
Intoxilyzer Serial # 80-000331 Location Vehicle #432 - Tr 2

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 104 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7sc

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
03/13/2009 11:42:14

Standard Lot#: 8112011
Last Changed By: C. CONLEY #34517

OAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:42:41
Cal Check	0.104	11:43:00
Air Blank	0.000	11:43:29

Sarah Cione

Operator's Signature