



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
9.31.2009

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 08/20/09 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 09/23/09 and 08/20/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


08/31/09

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 31, day of August, 2009.

A.M. Hernandez

Notary Public

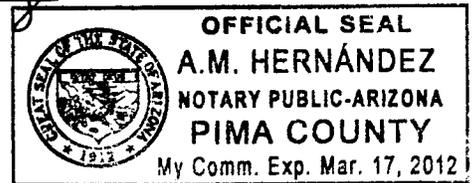

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

187
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R15-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

B4
8-31-09

QA SPECIALIST Cheryl Straub AGENCY Tucson Police Department

DATE 08/20/09 TIME 11:16

INTOXILYZER SERIAL # 80-000331 LOCATION Dut Squad HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

287
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
08/20/2009

11:18:04

TS4

Standard Lot#: 811201I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:19:02
Diagnostic Test	Pass	11:19:30
Air Blank	0.000	11:19:57
0.100 Cal Check	0.098	11:20:16
Air Blank	0.000	11:20:46
Subject Test	0.000	11:21:12
Air Blank	ABT*	11:21:20
Air Blank	0.000	11:21:48

*Sequence Aborted

Not a Successfully Completed Test Sequence

377
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
08/20/2009 11:22:05

BY

Standard Lot#: 811201I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	11:22:47
Air Blank	IPS**	11:22:56

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

487
C8

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
08/20/2009

11:23:14

BS4

Standard Lot#: 811201I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:24:11
Diagnostic Test	Pass	11:24:38
Air Blank	0.000	11:25:05
0.100 Cal Check	0.098	11:25:24
Air Blank	0.000	11:25:55
Subject Test	0.067*	11:29:03
Air Blank	0.000	11:29:35
Air Blank	0.000	11:30:02
Subject Test	XXX**	11:30:19
Air Blank	0.000	11:30:53

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

587
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000331

08/20/2009

11:31:17

Standard Lot#: 811201I

Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772

TUCSON PD

TS4

Test	g/210L	Time
Air Blank	0.000	11:31:44
Cal Check	0.098	11:32:03
Air Blank	0.000	11:32:33



Operator's Signature

687
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

TS4
8.31.09

QA Specialist Cheeryl STRAUB AGENCY Tucson Police Department
DATE 08/20/09 TIME 11:14
Intoxilyzer Serial # 80-000331 Location DUI Squad HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheeryl Straub

DPS Form Exh G-2 (Rev. 05-01)

787
CS


INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
08/20/2009 11:16:04

104

Standard Lot#: 811201I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:16:31
Cal Check	0.098	11:16:50
Air Blank	0.000	11:17:19

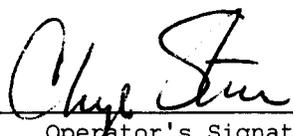

Operator's Signature