



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
5.2.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4/28/11 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 3/31/11 and 4/28/11  
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -

The instrument was returned to service on -

Sarah Cione

Sarah L. Cione

4/28/11


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 28, day of  
April, 20 11.

A.M. Hernandez

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

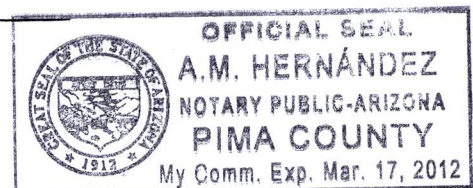

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7<sub>SC</sub>

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TSy  
5.2.2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 4/28/11 TIME 926

INTOXILYZER SERIAL # 80-000331 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES ✓ NO       

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000331  
04/28/2011

09:28:06

2 of 7<sub>sc</sub>

124

Standard Lot#: 654682  
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:28:49
Diagnostic Test	Pass	09:29:17
Air Blank	0.000	09:29:44
0.100 Cal Check	0.098	09:30:03
Air Blank	0.000	09:30:33
Subject Test	0.000	09:32:02
Air Blank	0.000	09:32:29
Five Minute Wait	ABT*	
Air Blank	0.000	09:33:10

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000331

04/28/2011

09:34:06

Standard Lot#: 654682

Last Changed By: C. CONLEY

#34517

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:34:50
Air Blank	IPS**	09:34:58

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000331  
04/28/2011

09:40:53

4 of 7<sub>sc</sub>

184

Standard Lot#: 654682  
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:41:41
Diagnostic Test	Pass	09:42:09
Air Blank	0.000	09:42:36
0.100 Cal Check	0.097	09:42:55
Air Blank	0.000	09:43:25
Subject Test	0.340*	09:46:33
Air Blank	0.000	09:47:05
Air Blank	0.000	09:47:33
Subject Test	XXX**	09:48:08
Air Blank	0.000	09:48:42

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**



INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000331

04/28/2011

09:49:33

Standard Lot#: 654682

Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:50:00
Cal Check	0.097	09:50:19
Air Blank	0.000	09:50:49

*Sarah Cione*  
Operator's Signature

5 of 7 SC

134

EXHIBIT G-2

6 of 7sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
5-2-2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 4/28/11 TIME 926  
Intoxilyzer Serial # 80-000331 Location DUT HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (Y) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000331

04/28/2011

09:26:32

-----  
Standard Lot#: 654682

Last Changed By: C. CONLEY #34517  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:26:59
Cal Check	0.097	09:27:18
Air Blank	0.000	09:27:47

*Sarah Cione*

Operator's Signature