



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

By
5.31.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/26/11 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 4/28/11 and 5/26/11
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -


Sarah Cione

Sarah L. Cione

5/30/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 30, day of

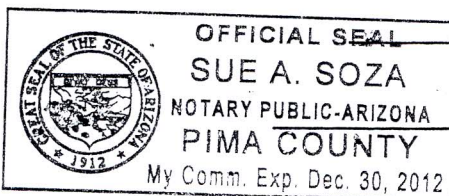

MAY, 20 11

Sue A. Soza

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

B6
5.31.2011

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 5/26/11 TIME 857
INTOXILYZER SERIAL # 80-000331 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 096 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7_{sc}

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000331

05/26/2011

08:59:49

BH

Standard Lot#: 654682

Last Changed By: C. CONLEY

#34517

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:00:33
Diagnostic Test	Pass	09:01:00
Air Blank	0.000	09:01:28
0.100 Cal Check	0.097	09:01:46
Air Blank	0.000	09:02:16
Subject Test	0.000	09:02:31
Air Blank	0.000	09:02:58
Five Minute Wait	ABT*	
Air Blank	0.000	09:03:35

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7_{sc}

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000331

05/26/2011

09:03:56

B4

Standard Lot#: 654682

Last Changed By: C. CONLEY

#34517

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:04:37
Air Blank	IPS**	09:04:45

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7_{sc}

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000331

05/26/2011

09:05:02

104

Standard Lot#: 654682

Last Changed By: C. CONLEY

#34517

QAS: S. CIONE

#50616

TUCSON PD

Operator: S. CIONE

#50616

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:05:59
Diagnostic Test	Pass	09:06:26
Air Blank	0.000	09:06:54
0.100 Cal Check	0.097	09:07:12
Air Blank	0.000	09:07:42
Subject Test	0.034*	09:10:50
Air Blank	0.000	09:11:22
Air Blank	0.000	09:11:49
Subject Test	XXX**	09:12:10
Air Blank	0.000	09:12:43

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000331

05/26/2011

09:17:24

Standard Lot#: 654682

Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:17:51
Cal Check	0.096	09:18:10
Air Blank	0.000	09:18:40

Sarah Cione

Operator's Signature

5 of 7_{sc}

BSH

EXHIBIT G-2

6 of 7_{sc}

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

BK
5.31.2011

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 5/26/11 TIME 857
Intoxilyzer Serial # 80-000331 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000331

05/26/2011

08:57:45

Standard Lot#: 654682

Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:58:12
Cal Check	0.097	08:58:31
Air Blank	0.000	08:59:00

Sarah Cione

Operator's Signature