



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
9.28.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9-15-11 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 8-18-11 and 9-15-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

9-23-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 23, day of September, 2011.

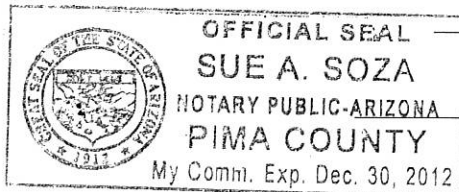

Sue A. Soza

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 9-15-11 TIME 0740

INTOXILYZER SERIAL # 80-000331 LOCATION ODW DL1 HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ✓ NO       

COMMENTS None.

SIGNATURE Clare Conley

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1509  
W

TS4  
9.28.2011

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 9-15-11 TIME 0740

INTOXILYZER SERIAL # 40-000331 LOCATION ODW Du1 HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.

Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.

Not a Successfully Completed Test Sequence printed

Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC

Result: 0. 098 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS Deficient sample was not performed

SIGNATURE Clare Conley



INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000331

09/15/2011

07:42:07

Standard Lot#: 04911100A1

Last Changed By: J. OLIVARES

#36552

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:42:58
Diagnostic Test	Pass	07:43:26
Air Blank	0.000	07:43:53
0.100 Cal Check	0.099	07:44:12
Air Blank	0.000	07:44:42
Subject Test	0.000	07:45:34
Air Blank	ABT*	07:45:40
Air Blank	0.000	07:46:09

\*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000331

09/15/2011

07:46:31

Standard Lot#: 04911100A1

Last Changed By: J. OLIVARES

#36552

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:47:02
Air Blank	IPS**	07:47:10

\*RFI Detect

\*\*Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000331

09/15/2011

07:47:49

Standard Lot#: 04911100A1

Last Changed By: J. OLIVARES

#36552

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:48:36
Diagnostic Test	Pass	07:49:04
Air Blank	0.000	07:49:31
0.100 Cal Check	0.099	07:49:50
Air Blank	0.000	07:50:20
Subject Test	XXX*	07:51:01
Air Blank	0.000	07:51:34

\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000331

09/15/2011

07:52:17

5 of 9  
cu

-----  
Standard Lot#: 04911100A1

Last Changed By: J. OLIVARES #36552  
-----

QAS: C. CONLEY #34517

TUCSON PD  
-----

TS4

Test	g/210L	Time
Air Blank	0.000	07:52:45
Cal Check	0.098	07:53:03
Air Blank	0.000	07:53:34

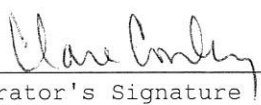
  
Operator's Signature



EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 9-15-11 TIME 0740  
Intoxilyzer Serial # 80-800331 Location ODW DUI HR

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 150 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000331

09/15/2011

07:40:41

7569  
CW

-----  
Standard Lot#: 04911100A1

Last Changed By: J. OLIVARES #36552  
-----

TS4

QAS: C. CONLEY #34517

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	07:41:09
Cal Check	0.099	07:41:27
Air Blank	0.000	07:41:57



Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

8569  
Cu

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

TS4  
9.28.2011

Date: 08-26-11  
QAS: OLIVARES  
Intoxilyzer Serial Number: 80-000331  
New Lot Number: 049 11100A1  
New Expiration Date: 03-01-13

- (L) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (4) 2. Remove the old gas tank and install the new one.
- (4) 3. Hit the "ESC" key twice and swipe the QAS card.
- (4) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
Enter Target Value: 0.100  
Enter Control Lot #: From the tank label  
Enter Exp Date: From the tank label (mm/dd/yyyy)  
Edit info Y/N? Y or N  
Update Std Info? Y ENTER  
"Please wait saving..."
- (7) 5. Hit the "ESC" key twice.
- (7) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
Select "D" for dry control test. After selection is made press ENTER.
- (8) 7. Air blank completed.
- (7) 8. Calibration check completed. Test results 0. 099 AC.
- (✓) 9. Air blank completed.
- (x) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE

[Signature] 36552

9569  
cm

TS4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-J01331  
06/27/2011 02:05:03

Standard Lot#: 04911100A1  
Last Changed By: J. OLIVARES #36552

OAS: J. OLIVARES #36552  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	02:05:31
Cal Check	0.099	02:05:49
Air Blank	0.000	02:06:19

  
Operator's Signature 36552