



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TSy
4-10-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/30/12 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 3/1/12 and 3/30/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


4/9/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 9th, day of April, 2012.

Elma Escarcega

Notary Public

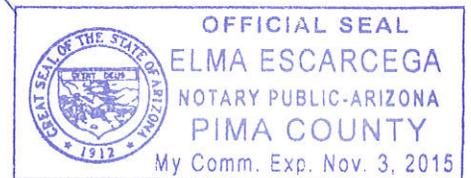

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 7 sc

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

By
4.10.2012

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 3/30/12 TIME 1356

INTOXILYZER SERIAL # 80-000331 LOCATION DUIHQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Sarah Cione

2 of 7sc

TSH

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
03/30/2012 13:57:57

Standard Lot#: 04911100A1
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:58:40
Diagnostic Test	Pass	13:59:08
Air Blank	0.000	13:59:35
0.100 Cal Check	0.099	13:59:54
Air Blank	0.000	14:00:23
Subject Test	0.000	14:00:43
Air Blank	0.000	14:01:10
Five Minute Wait	ABT*	
Air Blank	0.000	14:01:49

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7 sc
p4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
03/30/2012 14:01:55

Standard Lot#: 04911100A1
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:02:42
Diagnostic Test	Pass	14:03:09
Air Blank	0.000	14:03:37
0.100 Cal Check	RFI*	14:03:50
Air Blank	IPS**	14:03:58

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7_{sc}

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
03/30/2012 14:04:05

Standard Lot#: 04911100A1
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:04:50
Diagnostic Test	Pass	14:05:17
Air Blank	0.000	14:05:45
0.100 Cal Check	0.099	14:06:03
Air Blank	0.000	14:06:33
Subject Test	0.045*	14:09:41
Air Blank	0.000	14:10:13
Air Blank	0.000	14:10:40
Subject Test	XXX**	14:10:58
Air Blank	0.000	14:11:31

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 7_{sc}

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
03/30/2012 14:11:45

TS4

Standard Lot#: 04911100A1
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:12:13
Cal Check	0.098	14:12:31
Air Blank	0.000	14:13:02

Sarah Cione
Operator's Signature

6 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

TS4
4.10.2012

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 3/30/12 TIME 1356
Intoxilyzer Serial # 80-000331 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0.100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0.098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
03/30/2012 13:56:34

7 of 7 sc

Standard Lot#: 04911100A1
Last Changed By: J. OLIVARES #36552

TS4

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:57:01
Cal Check	0.098	13:57:20
Air Blank	0.000	13:57:50

Sarah Cione

Operator's Signature