



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

T54
5.7-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4/27/12 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 3/30/12 and 4/27/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


5/3/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 3rd day of May, 2012.

Elma Escarcega

Notary Public

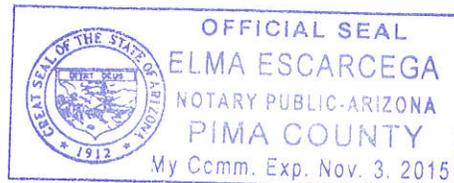

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 7 SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

TS4
5.7.2012

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 4/27/12 TIME 1015

INTOXILYZER SERIAL # 80-000331 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
04/27/2012

10:16:32

2 of 75c
B4

Standard Lot#: 04911100A1
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:17:15
Diagnostic Test	Pass	10:17:42
Air Blank	0.000	10:18:10
0.100 Cal Check	0.099	10:18:28
Air Blank	0.000	10:18:59
Subject Test	0.000	10:19:28
Air Blank	0.000	10:19:55
Five Minute Wait	ABT*	
Air Blank	0.000	10:20:29

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
04/27/2012

10:20:34

3 of 7 sc
BSH

Standard Lot#: 04911100A1
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:21:09
Air Blank	IPS**	10:21:16

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7 sc

By

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
04/27/2012 10:21:21

Standard Lot#: 04911100A1
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:22:02
Diagnostic Test	Pass	10:22:29
Air Blank	0.000	10:22:56
0.100 Cal Check	0.099	10:23:15
Air Blank	0.000	10:23:45
Subject Test	0.077*	10:26:54
Air Blank	0.000	10:27:25
Air Blank	0.000	10:27:53
Subject Test	XXX**	10:28:14
Air Blank	0.000	10:28:47

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 7 sc

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
04/27/2012 10:29:01

TS4

Standard Lot#: 04911100A1
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:29:28
Cal Check	0.098	10:29:47
Air Blank	0.000	10:30:17

Sarah Cione

Operator's Signature

EXHIBIT G-2

6 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

T34

5-7-2012

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 4/27/12 TIME 1015
Intoxilyzer Serial # 80-000331 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
04/27/2012 10:15:10


7 of 7 sc

TS4

Standard Lot#: 04911100A1
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:15:37
Cal Check	0.099	10:15:56
Air Blank	0.000	10:16:26

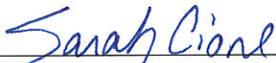

Operator's Signature