



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
9-20-2012

I hereby certify that these documents are a true and correct copy of the
periodic maintenance* for the Intoxilyzer 8000 Serial Number
80-000331 maintained by the Tucson Police Department
Crime Laboratory pursuant to the requirements of the Arizona Department
of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following
dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2,
for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/13/12 Cal 1 0.096 AC Cal 2 0.095 AC

Between the dates of 8/16/12 and 9/13/12
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on 9/13/12.
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione

9/17/12


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 17th, day of
September, 2012.

Elma Escarcega

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

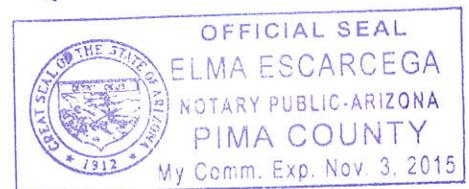

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the
information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 95c

ARIZONA DEPARTMENT OF PUBLIC SAFETY

TS4

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

9-20-2012

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 9/13/12 TIME 1012

INTOXILYZER SERIAL # 80-000331 LOCATION DUI HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 095 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 95c

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
09/13/2012 10:14:16

Standard Lot#: 04911100A1
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:15:00
Diagnostic Test	Pass	10:15:27
Air Blank	0.000	10:15:54
0.100 Cal Check	0.097	10:16:13
Air Blank	0.000	10:16:43
Subject Test	0.000	10:17:36
Air Blank	0.000	10:18:04
Five Minute Wait	ABT*	
Air Blank	0.000	10:18:41

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 9 sc

T34

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
09/13/2012 10:22:01

Standard Lot#: 04911100A1
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:22:34
Air Blank	IPS**	10:22:44

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 9 sc

T24

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
09/13/2012 10:27:21

Standard Lot#: 04911100A1
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:28:25
Diagnostic Test	Pass	10:28:53
Air Blank	0.000	10:29:20
0.100 Cal Check	0.097	10:29:39
Air Blank	0.000	10:30:09
Subject Test	0.103*	10:33:17
Air Blank	0.000	10:33:49
Air Blank	0.000	10:34:17
Subject Test	XXX**	10:34:49
Air Blank	0.000	10:35:22

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 9 sc


TSh

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
09/13/2012 10:36:01

Standard Lot#: 04911100A1
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:36:28
Cal Check	0.095	10:36:47
Air Blank	0.000	10:37:17

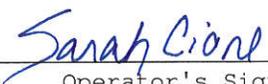

Operator's Signature

6 of 95c

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

TSG
9.20.2012

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 9/13/12 TIME 1012
Intoxilyzer Serial # 80-000331 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
09/13/2012 10:12:38


7 of 9 JC

TS4

Standard Lot#: 04911100A1
Last Changed By: J. OLIVARES #36552

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:13:05
Cal Check	0.096	10:13:24
Air Blank	0.000	10:13:54

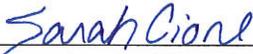

Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

8 of 9 sc

INTOXILYZER SERVICE RECORD

TS4
9.20.2012

Date: 9/13/12

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-000331

Reason For Service: Calibration checks starting to run -5%, Preventative maintenance.

Date Out of Service: 9/13/12

Calibration Check completed prior to removal from service

Date: 9/13/12

SQAP completed prior to removal from service

Date: 9/13/12

Corrective Action: Sent to CMI, Inc.

9 of 9 sc

TS4



Service Evaluation Form

This form MUST be completed and enclosed with instrument to be serviced.
Failure to complete and return this form may cause delays in service.

(Note: please ship items in their original shipping container(s) or a similar protective box.)

1. Contact information: Customer Number 105075 (contact Customer Service)
 Name Terry Gallegos Phone: (520) 837-7677
 Fax: (520) 791-4642 Email: Terry.Gallegos@tucsonaz.gov

2. Bill to Address: Ship to Address:
Finance / Accounts Payable Tucson Police Department Crime Lab
Po Box 27450 1306 W. Miracle Mile
Tucson, AZ 85726-7450 Tucson, AZ 85705

3. Serial Number: 80-000331 Instrument Model: Intox 8000

4. Detailed Description of Problem:
Preventative maintenance. Calibration checks - 5%.

Hazardous Material Warning! - DO NOT return gas cylinders with instrument!

I Authorize Repairs Up To: All \$250 \$500 \$750 Other \$ _____
 Purchase Order Number (attach a copy of P.O. if applicable) 34770

Authorized By:

Terry Gallegos Toxicology Coordinator
 Name (Please Print) Title
[Signature] 9.20.2012
 Signature Date

Ship to:
 CMI, Inc.
 Attn: Service Dept.
 316 East Ninth Street
 Owensboro, KY 42303

No, please send estimate before repairs are made.
 Note: An estimate will be faxed before performing any repairs and may cause delays in service.
 An evaluation fee (\$79.00 or actual costs) will apply to estimates that are not repaired.