



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
10-31-2013

I hereby certify that these documents are a true and correct copy of the
periodic maintenance* for the Intoxilyzer 8000 Serial Number
80-000331 maintained by the Tucson Police Department
Crime Laboratory pursuant to the requirements of the Arizona Department
of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following
dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2,
for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/10/13 Cal 1 0.097 AC Cal 2 0.095 AC

Between the dates of 9/12/13 and 10/10/13
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


10/21/13

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 21st, day of
October, 2013.

Edma Escarcega

Notary Public

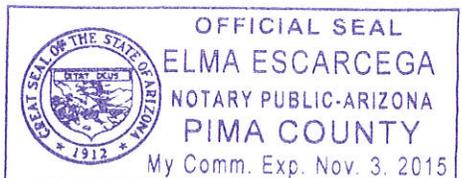

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7 sc
TS4
10-31-2013

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 10/10/13 TIME 958

INTOXILYZER SERIAL # 80-000331 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 095 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000331
Core Version: 8105.48
10/10/2013

10:01:04

2 of 7_{sc}
TS4

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:01:54
Diagnostic Test	Pass	10:02:28
Air Blank	0.000	10:02:56
0.100 Cal Check	0.096	10:03:17
Air Blank	0.000	10:03:46
Subject Test	0.000	10:04:13
Air Blank	0.000	10:04:41
Wait	ABT*	
Air Blank	0.000	10:05:17

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000331
Core Version: 8105.48
10/10/2013

10:08:58

3 of 7 sc
TS4

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:09:57
Diagnostic Test	Pass	10:10:31
Air Blank	0.000	10:10:59
0.100 Cal Check	0.096	10:11:19
Air Blank	0.000	10:11:49
Subject Test	RFI*	10:12:24
Air Blank	IPS**	10:12:31

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000331
Core Version: 8105.48
10/10/2013

10:12:37

4 of 7 sc
T84

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:13:36
Diagnostic Test	Pass	10:14:10
Air Blank	0.000	10:14:37
0.100 Cal Check	0.097	10:14:58
Air Blank	0.000	10:15:28
Subject Test	0.215*	10:18:35
Air Blank	0.000	10:19:08
Air Blank	0.000	10:19:35
Subject Test	XXX**	10:20:13
Air Blank	0.000	10:20:47

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000331
Core Version: 8105.48
10/10/2013 10:21:13

5 of 7
SC
TSK

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:21:42
0.100 Cal Check	0.095	10:22:02
Air Blank	0.000	10:22:32

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

6 of 7
sc
TS4
10-31-2013

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 10/10/13 TIME 958
Intoxilyzer Serial # 80-000331 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON PD

Serial Number: 80-000331

Core Version: 8105.48

10/10/2013

09:58:12

7 of 7 SC
TS4

Standard Lot#: 30412100A1

Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616

TUCSON PD

Operator: S. CIONE #50616

TUCSON PD

Test g/210L Time

Air Blank 0.000 09:58:40
0.100 Cal Check 0.097 09:59:01
Air Blank 0.000 09:59:31

Preventative Maintenance
Performed -- Timer Reset