



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
2.12.2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/30/14 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 1/3/14 and 1/30/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


2/11/14

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 11th, day of February, 2014.

Elma Escarcega

Notary Public

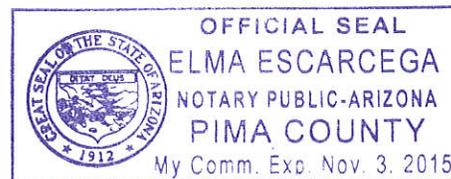

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
sc
By
2-12-2014

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 1/30/14 TIME 837

INTOXILYZER SERIAL # 80-000331 LOCATION DUI HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7
sc
TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000331
Core Version: 8105.48
01/30/2014 08:39:34

Standard Lot#: 2591310092
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:40:25
Diagnostic Test	Pass	08:41:01
Air Blank	0.000	08:41:29
0.100 Cal Check	0.101	08:41:49
Air Blank	0.000	08:42:19
Subject Test	0.000	08:42:52
Air Blank	0.000	08:43:19
Wait	ABT*	
Air Blank	0.000	08:44:01

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000331
Core Version: 8105.48
01/30/2014 08:44:25

3 of 7
sc
TSh

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:45:24
Diagnostic Test	Pass	08:46:01
Air Blank	0.000	08:46:28
0.100 Cal Check	0.101	08:46:49
Air Blank	0.000	08:47:19
Subject Test	RFI*	08:47:53
Air Blank	IPS**	08:48:05

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

INTOXALYZER 8000
Location: TUCSON PD
Serial Number: 80-000331
Core Version: 8105.48
01/30/2014 08:52:17

4 of 7 sc
TS4

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:53:15
Diagnostic Test	Pass	08:53:51
Air Blank	0.000	08:54:19
0.100 Cal Check	0.099	08:54:39
Air Blank	0.000	08:55:09
Subject Test	0.133*	08:58:17
Air Blank	0.000	08:58:50
Air Blank	0.000	08:59:17
Subject Test	XXX**	08:59:55
Air Blank	0.000	09:00:30

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully
Completed
Test Sequence

5 of 7_{sc}
TSH

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000331
Core Version: 8105.48
01/30/2014 09:04:37

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:05:05
0.100 Cal Check	0.099	09:05:26
Air Blank	0.000	09:05:56
Air Blank	0.000	09:08:24

*Input Timed Out

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7_{sc}
TSH
2.12.2014

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 1/30/14 TIME 837
Intoxilyzer Serial # 80-000331 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7 sc
TS4

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000331
Core Version: 8105.48
01/30/2014 08:37:33

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

OAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:38:01
0.100 Cal Check	0.100	08:38:22
Air Blank	0.000	08:38:52

Preventative Maintenance
Performed -- Timer Reset
